

New York State Police Crime Laboratory System
BREATH TESTING

Technical & Administrative Review Checklist - Simulator Solution

Class: Form
ID: Breath Form 4b (Version 5)

Effective Date: January 11, 2017
Page: 1 of 1

Solution Lot Number: 17310

Review Start Date: 10/25/17

Reviewer (initials): CA

Documentation:

- ☒ Handwritten initials on all pages
- ☒ Page number on all pages, total pages noted
- ☒ Solution lot number and bottle number on all solution pages
- ☒ Proper sample loading documented
- ☒ All pages of the technical record were reviewed for corrections; if applicable, corrections were documented appropriately
- ☒ Dates of analyses documented
- ☒ If abbreviations/symbols specific to the laboratory are used, they have been defined in section procedures
- NA Any deviations/departures from acceptability criteria have documented approval

Quality Control Data:

- ☒ Correct version of the test method is indicated on the QC Worksheet
- ☒ Required control samples analyzed, expiration dates not exceeded
- ☒ Acceptance criteria met
- ☒ QC results are accurately entered in Tox Form 11 QC Spreadsheet

Technical Data:

- ☒ The technical records demonstrate compliance with test methods
- ☒ Data transfers are accurate
- ☒ Retention times are recorded and acceptable
- ☒ Results meet acceptability criteria

Simulator Solution Record (electronic):

- ☒ Solution specific information is accurate (e.g. lot number, expiration date, etc)
- ☒ Results are accurately reported
- ☒ Electronic signature applied
- ☒ File name accurate

Review Completed:

- ☒ Upon review of the technical records available, I independently reached the same conclusion(s) reported.

Reviewer (initial/date):

CA 10/25/17

Page 1 of 1

**New York State Police Crime Laboratory System
BREATH TESTING**

Simulator Solution Verification Results

Class: Form
ID: Breath Form 2 (Version 3)

Effective Date: June 1, 2015

Page: 1 of 1

Simulator Solution Level: 0.10 %

Target Concentration: 0.121 %

Acceptable Range: 0.117-0.125 %

Manufacturer: Guth Laboratories, Inc.

Analyst: Carol S Miller

Lot Number: 17310

Analysis Date: 10/5/2017

Bottle Number	% by wt	AVG
6	0.1184	0.1187
	0.1189	
	0.1188	
7	0.1188	0.1188
	0.1188	
	0.1187	
705	0.1193	0.1191
	0.1190	
	0.1191	
706	0.1191	0.1191
	0.1190	
	0.1191	
1390	0.1191	0.1191
	0.1190	
	0.1191	
1391	0.1194	0.1192
	0.1195	
	0.1188	
Grand Average		0.1190

pg 1 of 57 *cm*

**New York State Police Crime Laboratory System
BREATH TESTING**

Simulator Solution Verification QC Worksheet

Class: Form
ID: Breath Form 8 (Version 2)

Effective Date: **January 10, 2017**
Page: 1 of 1

Analysis Date: 10-05-17
Analyst Initials: Cam

Breath Instr 2 Version 5

Simulator Solution Lot# 17310

Expiration Date 03-25-18
Received Date 10-02-17
Level 0.10%
Bottle Number(s) 67, 705, 706, 1390, 1391

<u>Internal Standard</u>	<u>Lot #</u>	<u>Expiration Date</u>
n-propanol	<u>SWL092017</u>	<u>12-20-17</u>

<u>Controls</u>	<u>Lot #</u>	<u>Expiration Date</u>
0.08 % by wt	<u>FN10281510</u>	<u>11/2020</u>
0.10 % by wt	<u>FN08101601</u>	<u>08/2021</u>
0.20 % by wt	<u>FN03301601</u>	<u>04/2021</u>
0.30 % by wt	<u>NA</u>	<u>NA</u>
0.40 % by wt	<u>NA</u>	<u>NA</u>
0.50 % by wt	<u>NA</u>	<u>NA</u>
Composite	<u>SWL022117</u>	<u>02-21-18</u>

<u>Instrumentation</u>	<u>Instrument ID/Serial #</u>	<u>Calibration Exp Date</u>
Headspace GC	<u>H-TX-HS2</u>	<u>12-20-17</u>
Dilutor	<u>MD965E1022</u>	<u>03-31-18</u>

Test Method Analyst Review

Batch Accepted ☒ Batch Rejected ☐

Rejection reason/course of action: _____

Deviations/departures: _____

Comments: _____

Initial/Date: Cam 10-25-17

2 Cam

Approved By: Top Management

Print Date: 10/5/2017

TotalChrom Sequence File C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

Printed by : manager on: 10/5/2017 1:30:14 PM
 Created by : toxicology on: 8/6/2009 1:12:05 PM
 Edited by : manager on: 10/5/2017 1:30:12 PM
 Number of Times Edited : 50

Description:

Sequence File Header Information:

Number of Rows : 27
 Instrument Type : PE AutoSystem GC
 Injection Type : SINGLE
 Raw tokens channel A :
 Result tokens channel A :
 Modified tokens channel A :
 Raw tokens channel B :
 Result tokens channel B :
 Modified tokens channel B :

Sequence Sample Descriptions - Channel A

Row	Type	Name
✓1	Sample	COMPOSITE ✓
✓2	Sample	0.08% CONTROL ✓
✓3	Sample	0.10% CONTROL ✓
✓4	Sample	NEG CONTROL ✓
✓5	Sample	GUTH 17310 #6 ✓
✓6	Sample	GUTH 17310 #6 ✓
✓7	Sample	GUTH 17310 #6 ✓
✓8	Sample	GUTH 17310 #7 ✓
✓9	Sample	GUTH 17310 #7 ✓
✓10	Sample	GUTH 17310 #7 ✓
✓11	Sample	GUTH 17310 #705 ✓
✓12	Sample	GUTH 17310 #705 ✓
✓13	Sample	GUTH 17310 #705 ✓
✓14	Sample	0.10% CONTROL ✓
✓15	Sample	NEG CONTROL ✓
✓16	Sample	GUTH 17310 #706 ✓
✓17	Sample	GUTH 17310 #706 ✓
✓18	Sample	GUTH 17310 #706 ✓
✓19	Sample	GUTH 17310 #1390 ✓
✓20	Sample	GUTH 17310 #1390 ✓
✓21	Sample	GUTH 17310 #1390 ✓
✓22	Sample	GUTH 17310 #1391 ✓
✓23	Sample	GUTH 17310 #1391 ✓
✓24	Sample	GUTH 17310 #1391 ✓
✓25	Sample	0.10% CONTROL ✓
✓26	Sample	0.20% CONTROL ✓
✓27	Sample	NEG CONTROL ✓

10-05-17


KMM
10/5/17

3 CMM

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33274
Sample Name : COMPOSITE
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 1

Date : 10/5/2017 2:43:15 PM
Data Acquisition Time : 10/5/2017 2:38:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-001.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

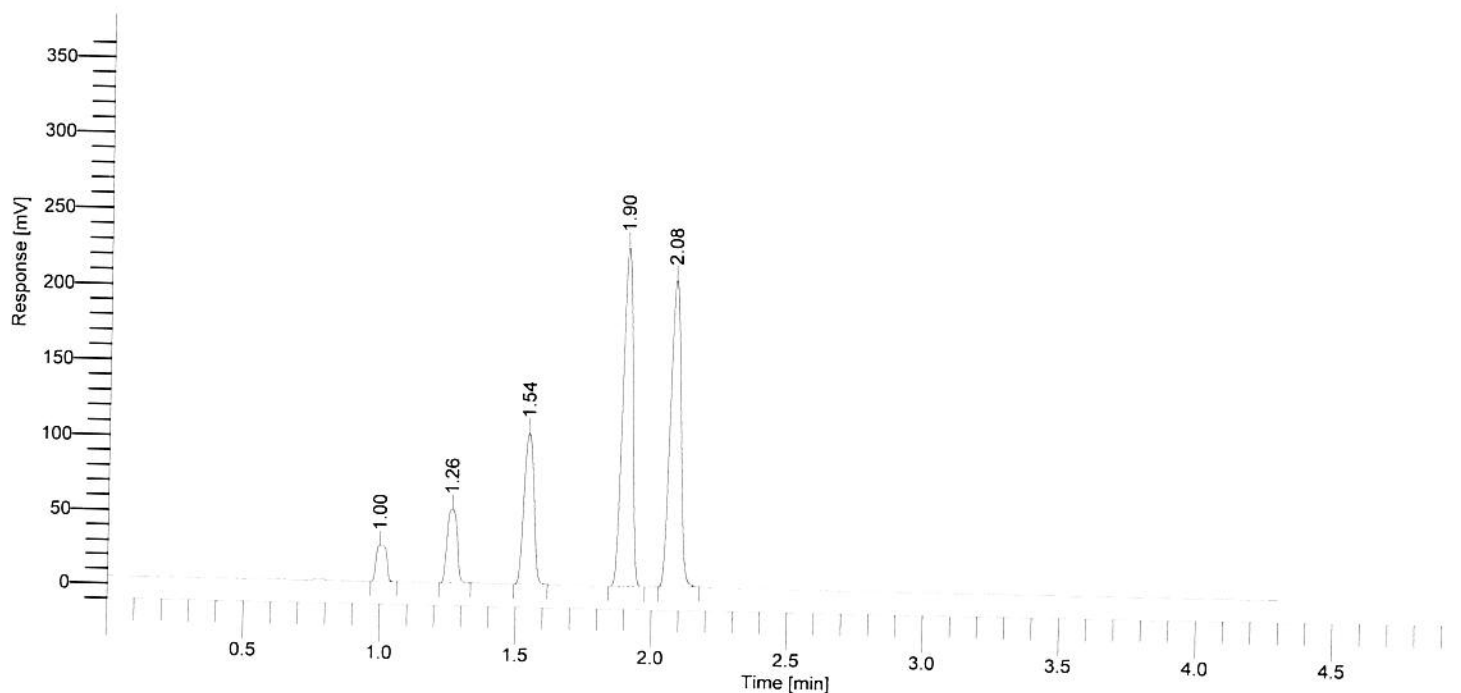

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0782	1.00	64443.2
ETHANOL	0.0774	1.26	125592.7
ISOPROPANOL	0.0777	1.54	255077.8
ACETONE	0.0761	1.90	572910.4
N-PROPANOL	-----	2.08	546374.0

400m

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33275
Sample Name : COMPOSITE
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 1

Date : 10/5/2017 2:43:16 PM
Data Acquisition Time : 10/5/2017 2:38:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-001.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

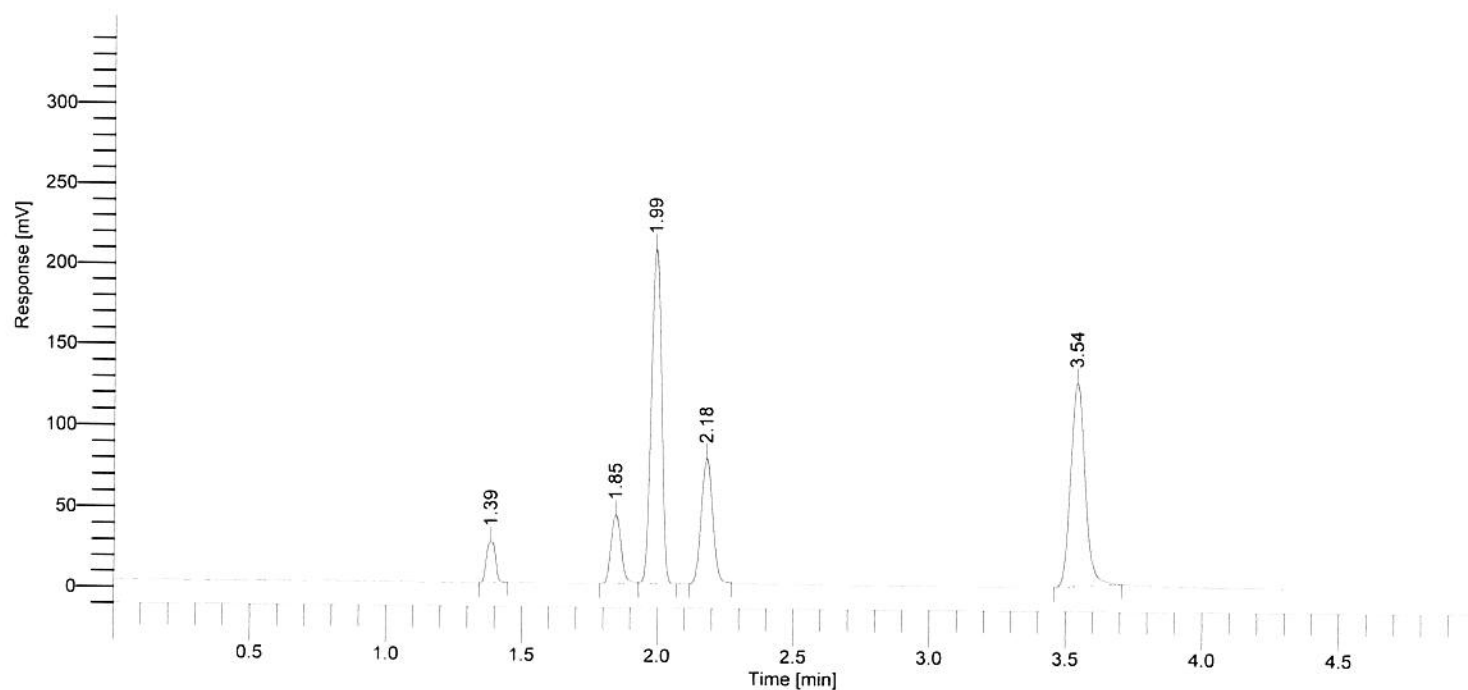

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
METHANOL	1.39
ETHANOL	1.85
ACETONE	1.99
ISOPROPANOL	2.18
N-PROPANOL	3.54

56m

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33276
Sample Name : 0.08% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 2

Date : 10/5/2017 2:48:13 PM
Data Acquisition Time : 10/5/2017 2:43:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-002.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

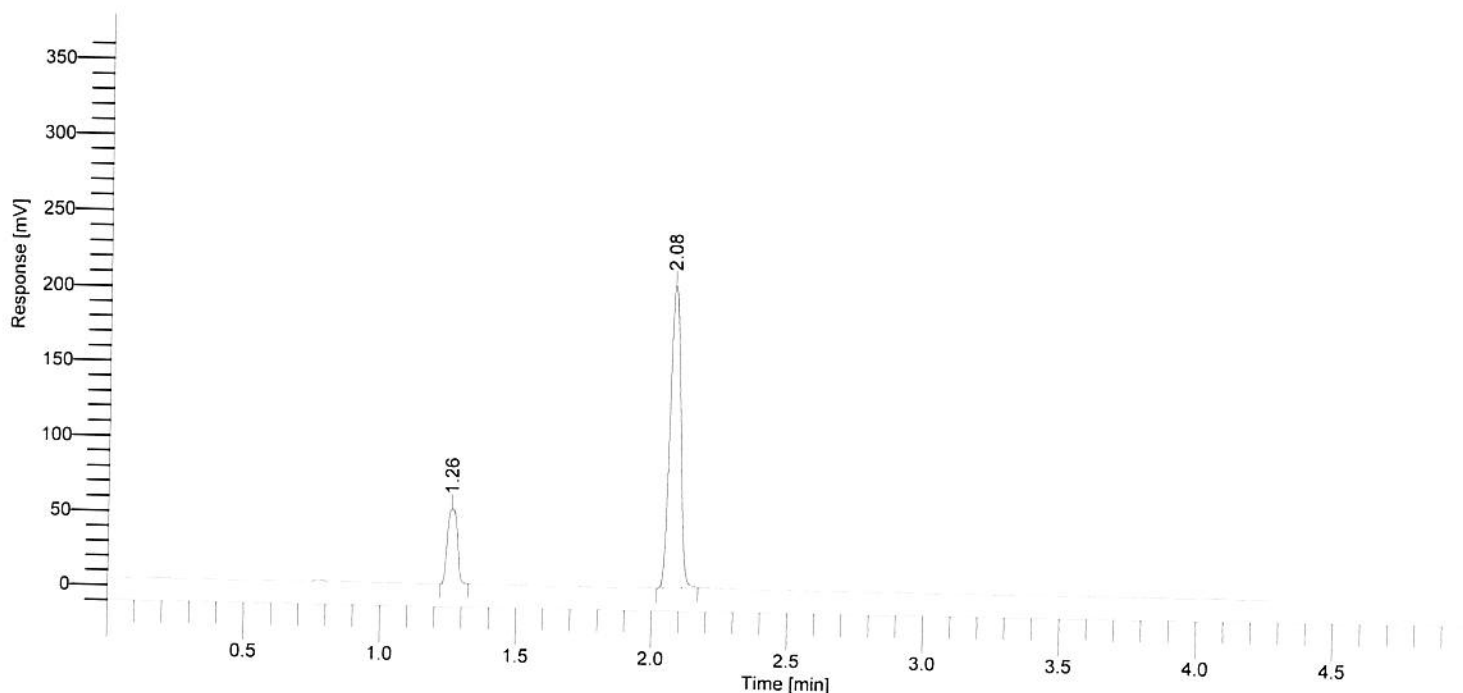

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.0781	1.26	125953.7
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	543243.8

6.2

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33277
Sample Name : 0.08% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 2

Date : 10/5/2017 2:48:15 PM
Data Acquisition Time : 10/5/2017 2:43:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-002.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

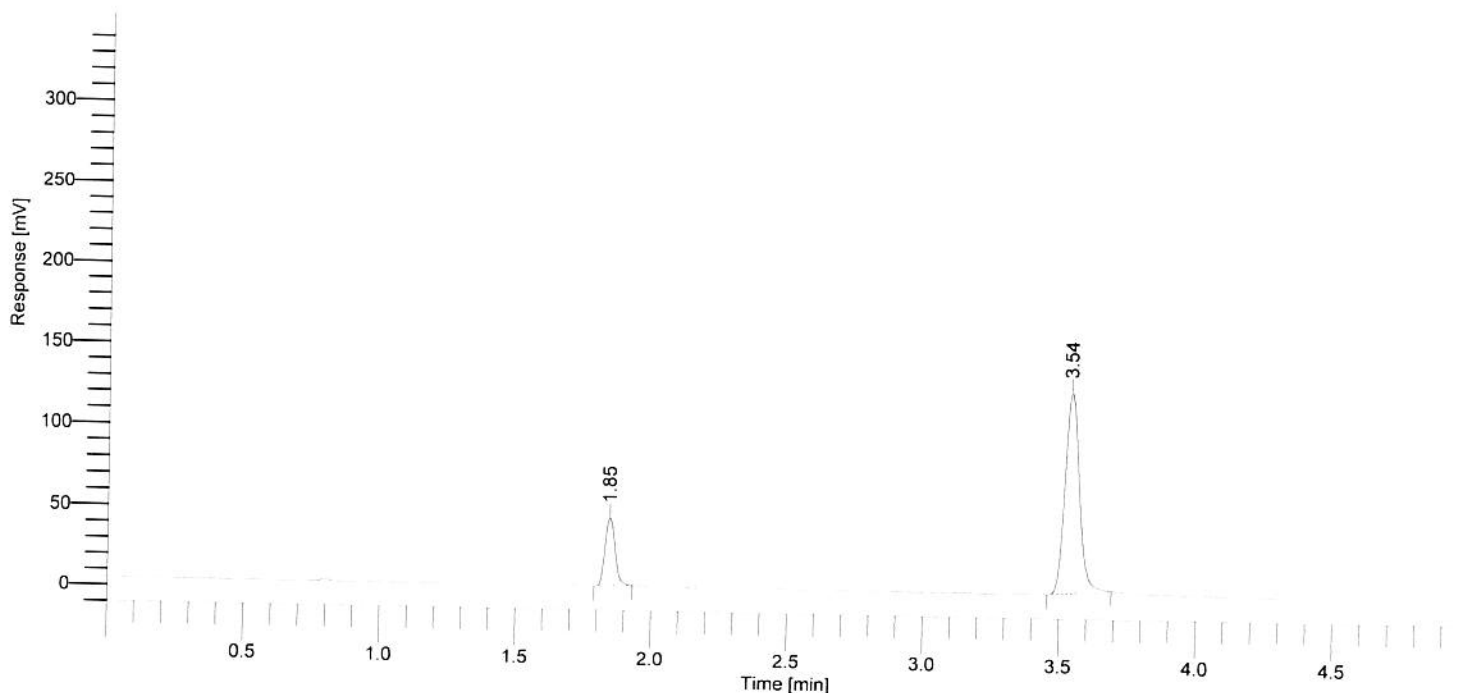

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

70m

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33278
Sample Name : 0.10% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 3

Date : 10/5/2017 2:53:13 PM
Data Acquisition Time : 10/5/2017 2:48:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-003.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

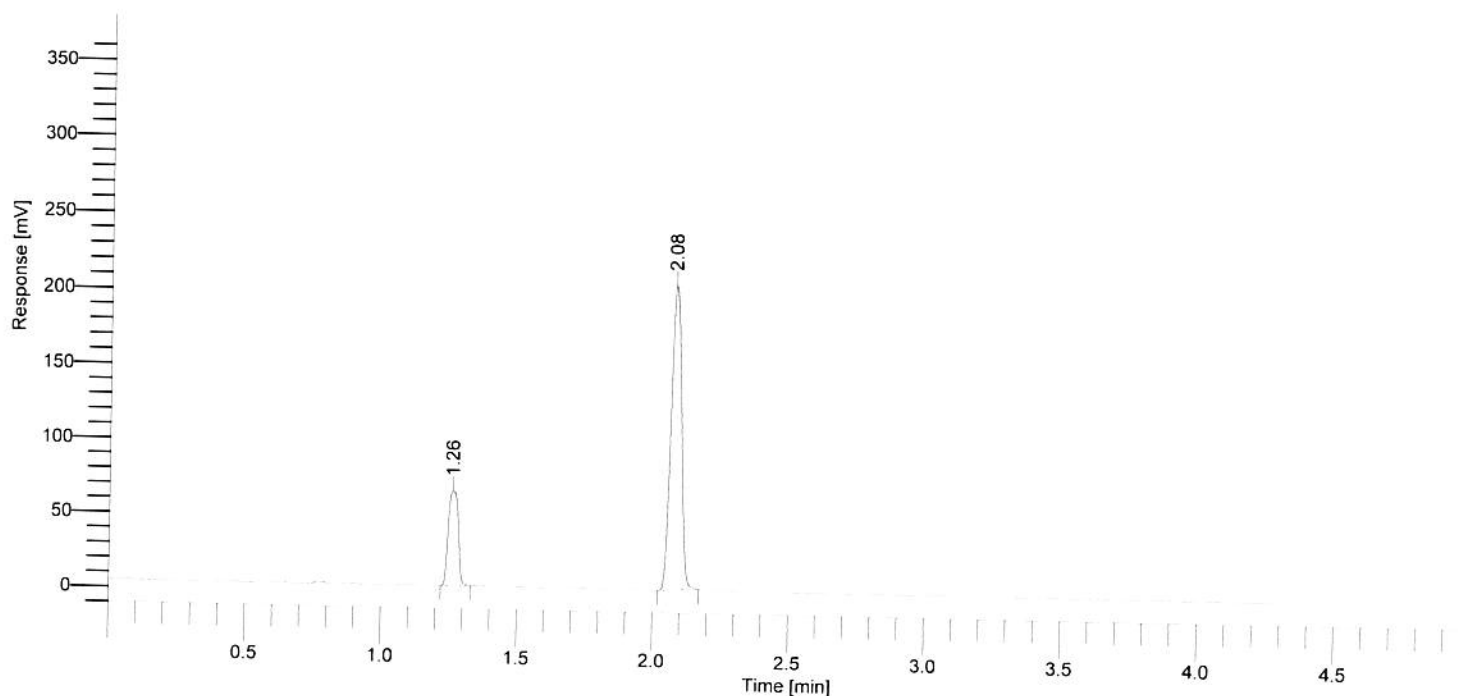

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.0979	1.26	158009.0
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	540539.9

8 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33279
Sample Name : 0.10% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 3

Date : 10/5/2017 2:53:15 PM
Data Acquisition Time : 10/5/2017 2:48:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-003.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

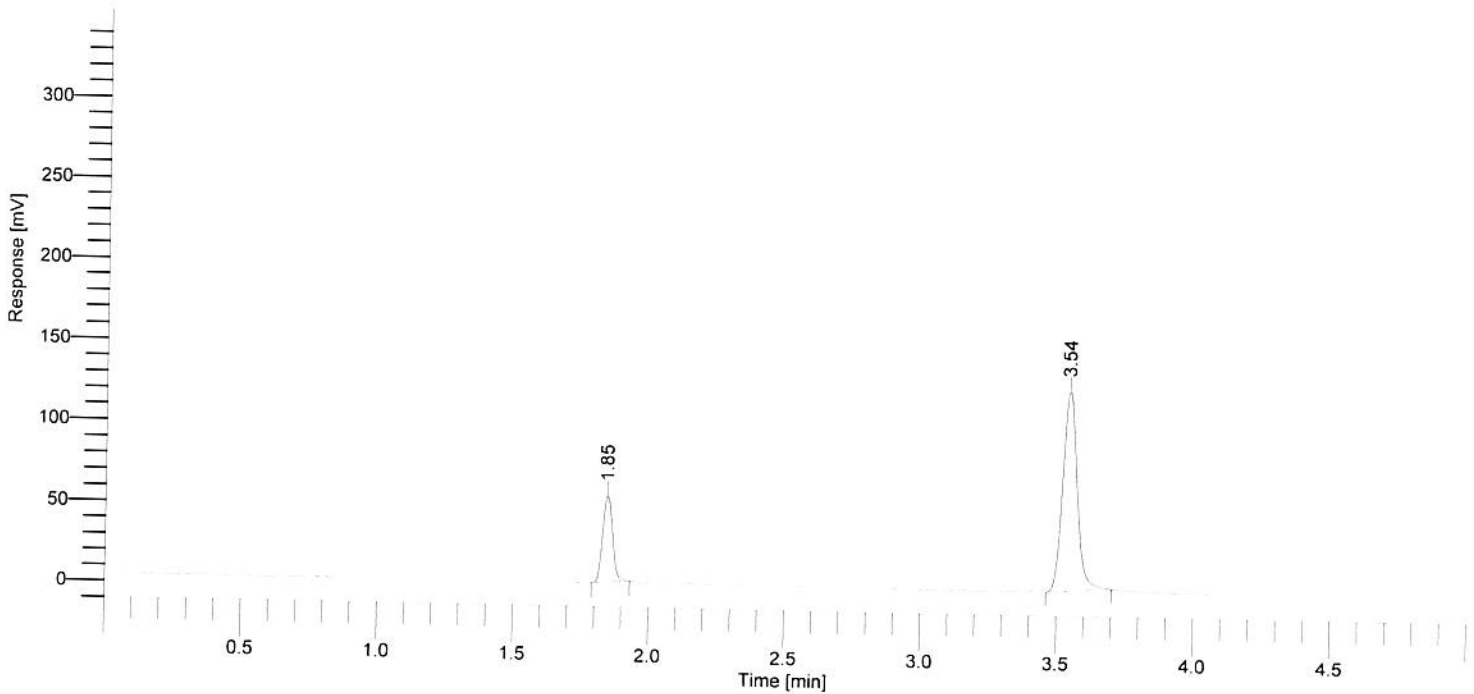

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

9 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33280
Sample Name : NEG CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 4

Date : 10/5/2017 2:58:13 PM
Data Acquisition Time : 10/5/2017 2:53:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-004.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

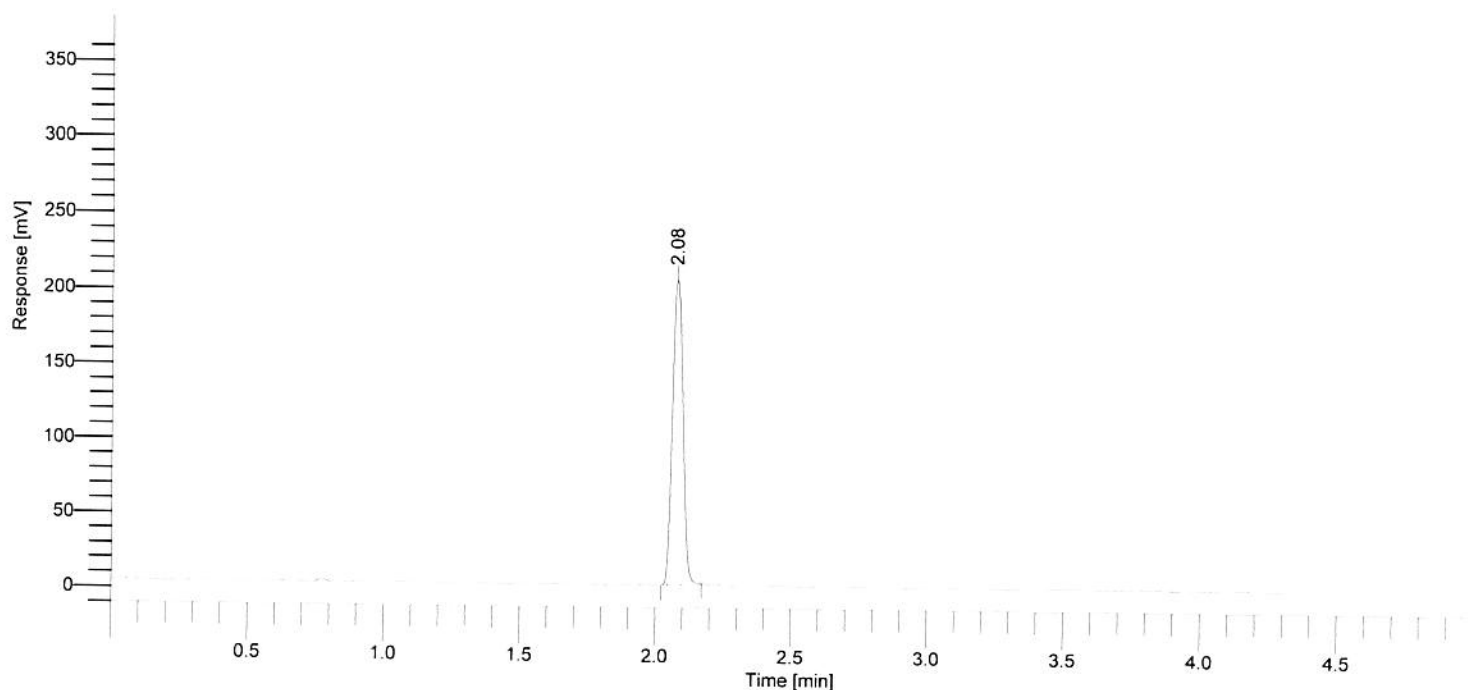

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu V \cdot s$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.0000	1.26	0.0
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	546033.4

10 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33281
Sample Name : NEG CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 4

Date : 10/5/2017 2:58:15 PM
Data Acquisition Time : 10/5/2017 2:53:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-004.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

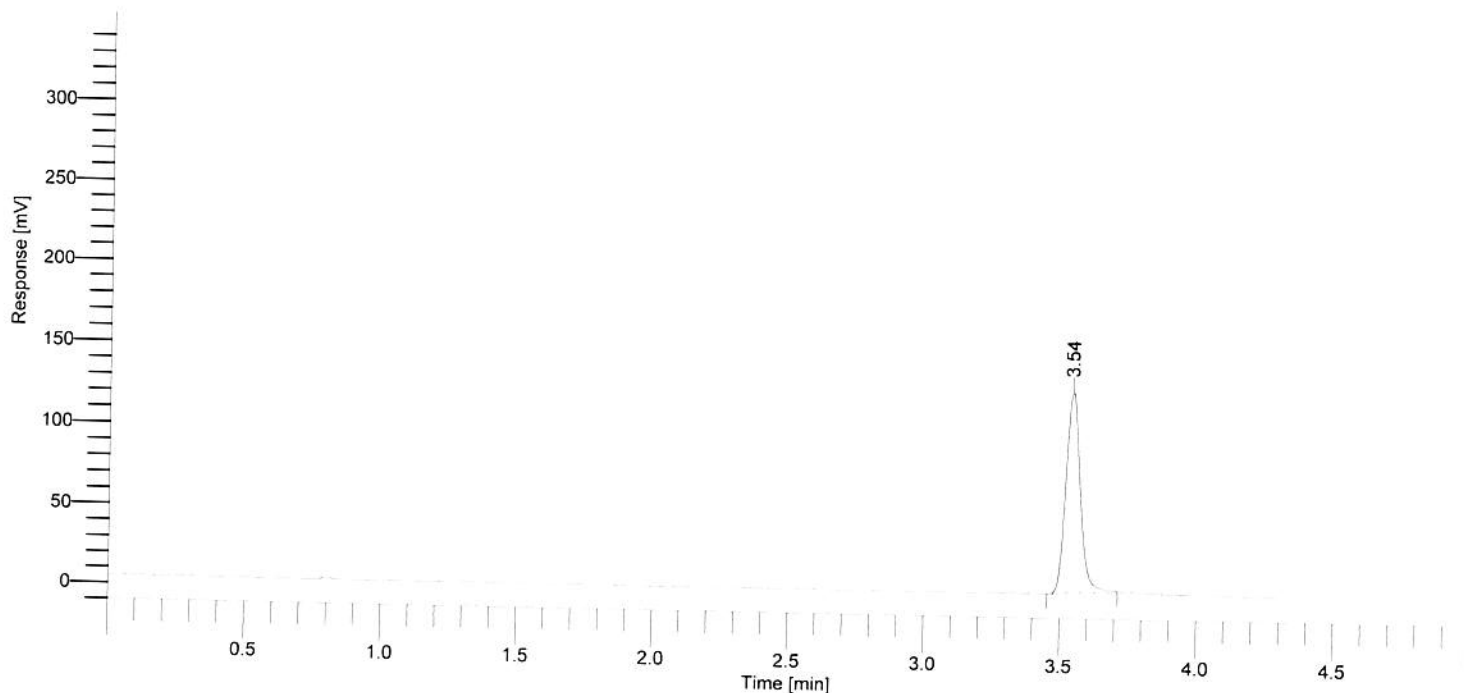

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
N-PROPANOL	3.54

11/6/17

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33282
Sample Name : GUTH 17310 #6
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 5

Date : 10/5/2017 3:03:13 PM
Data Acquisition Time : 10/5/2017 2:58:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-005.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

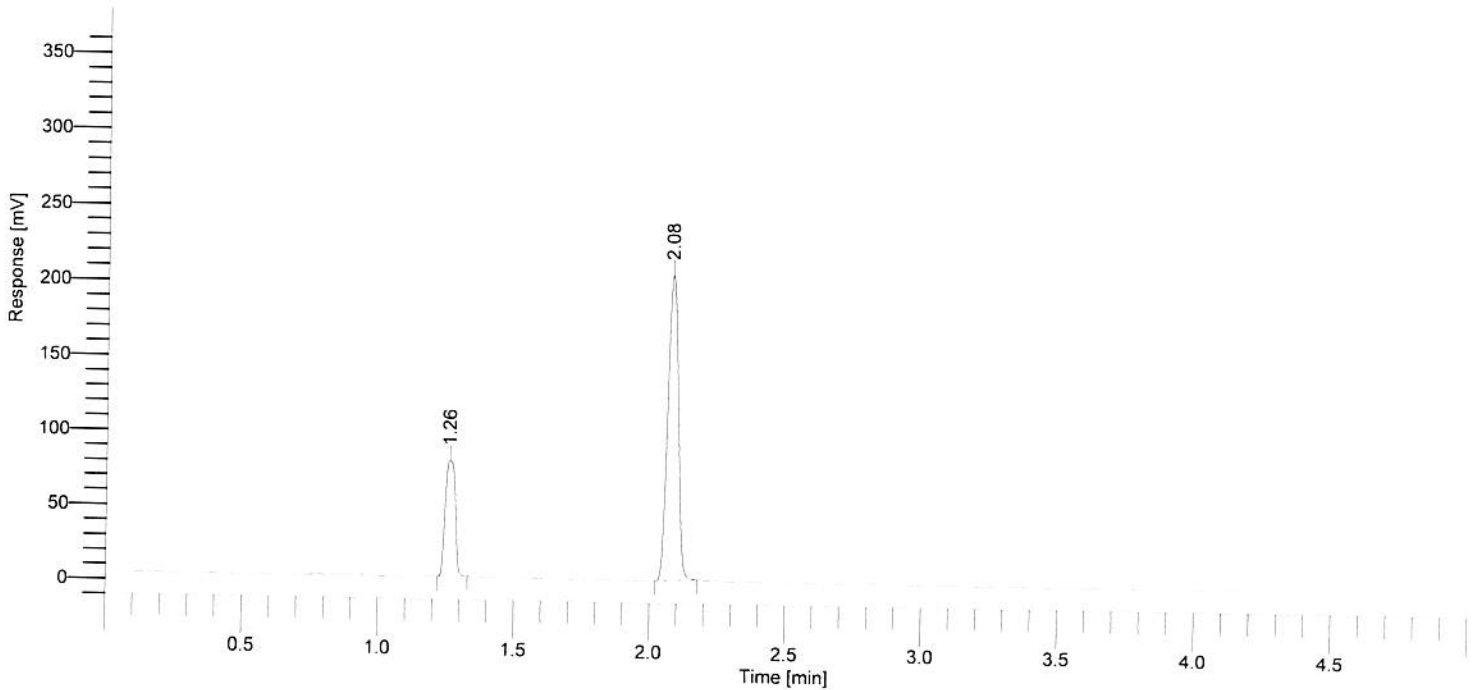

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1184	1.26	194628.2
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	548474.9

12 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33283
Sample Name : GUTH 17310 #6
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 5

Date : 10/5/2017 3:03:15 PM
Data Acquisition Time : 10/5/2017 2:58:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-005.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

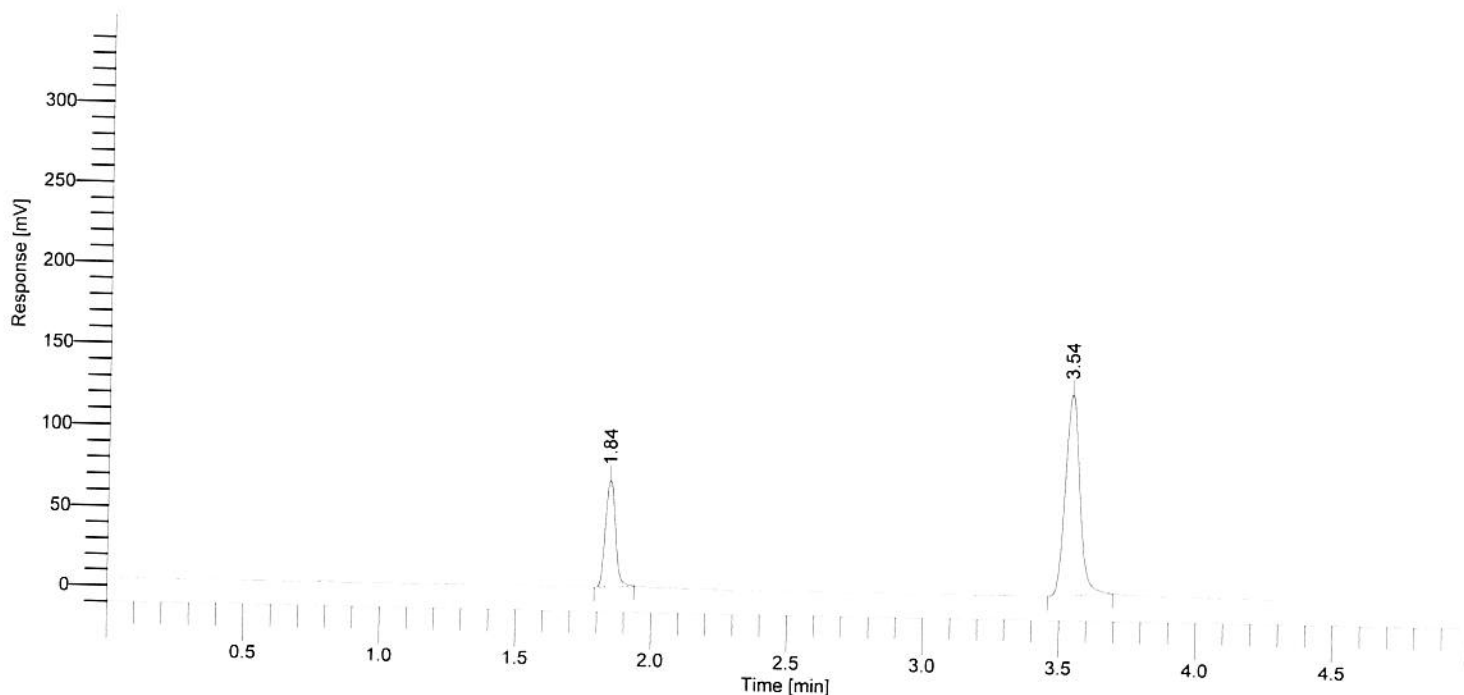

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.84
N-PROPANOL	3.54

13 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33284
Sample Name : GUTH 17310 #6
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 6

Date : 10/5/2017 3:08:13 PM
Data Acquisition Time : 10/5/2017 3:03:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-006.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

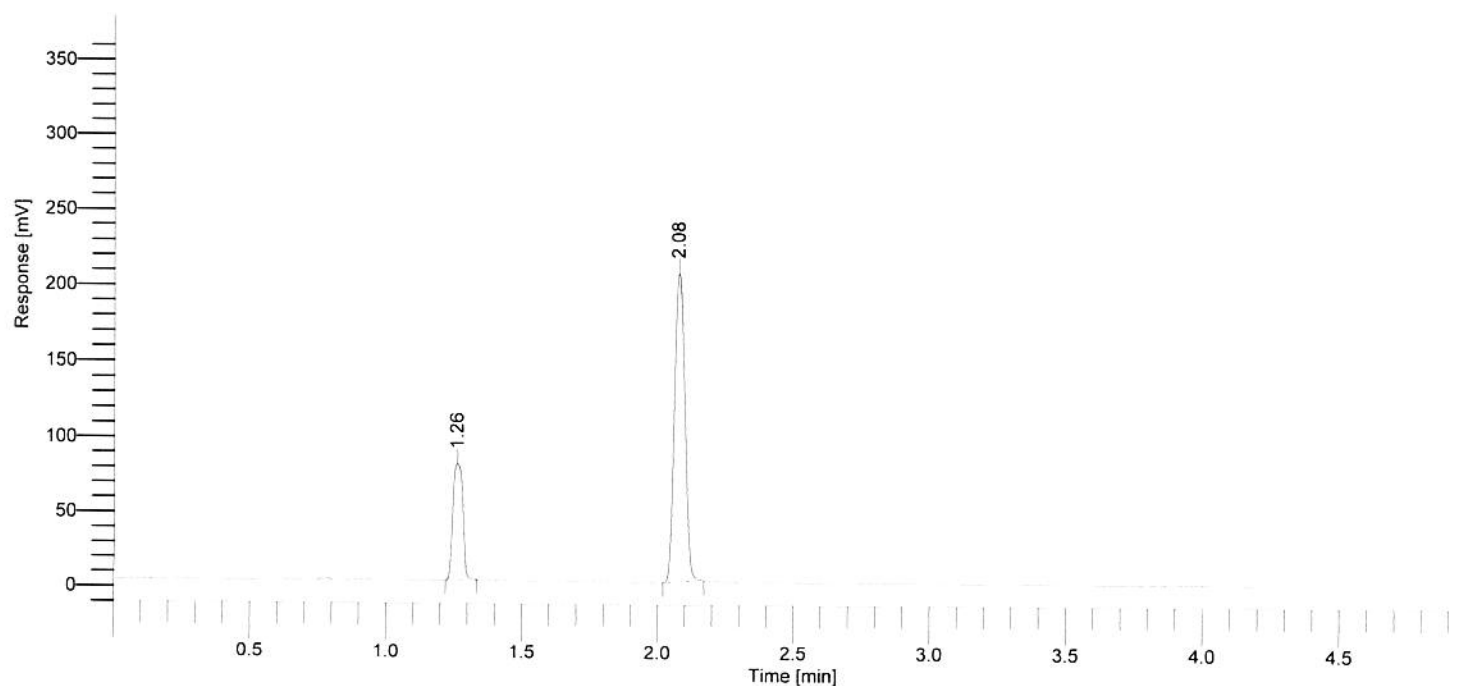

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1189	1.26	194807.2
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	546886.6

14 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33285
Sample Name : GUTH 17310 #6
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 6

Date : 10/5/2017 3:08:15 PM
Data Acquisition Time : 10/5/2017 3:03:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-006.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

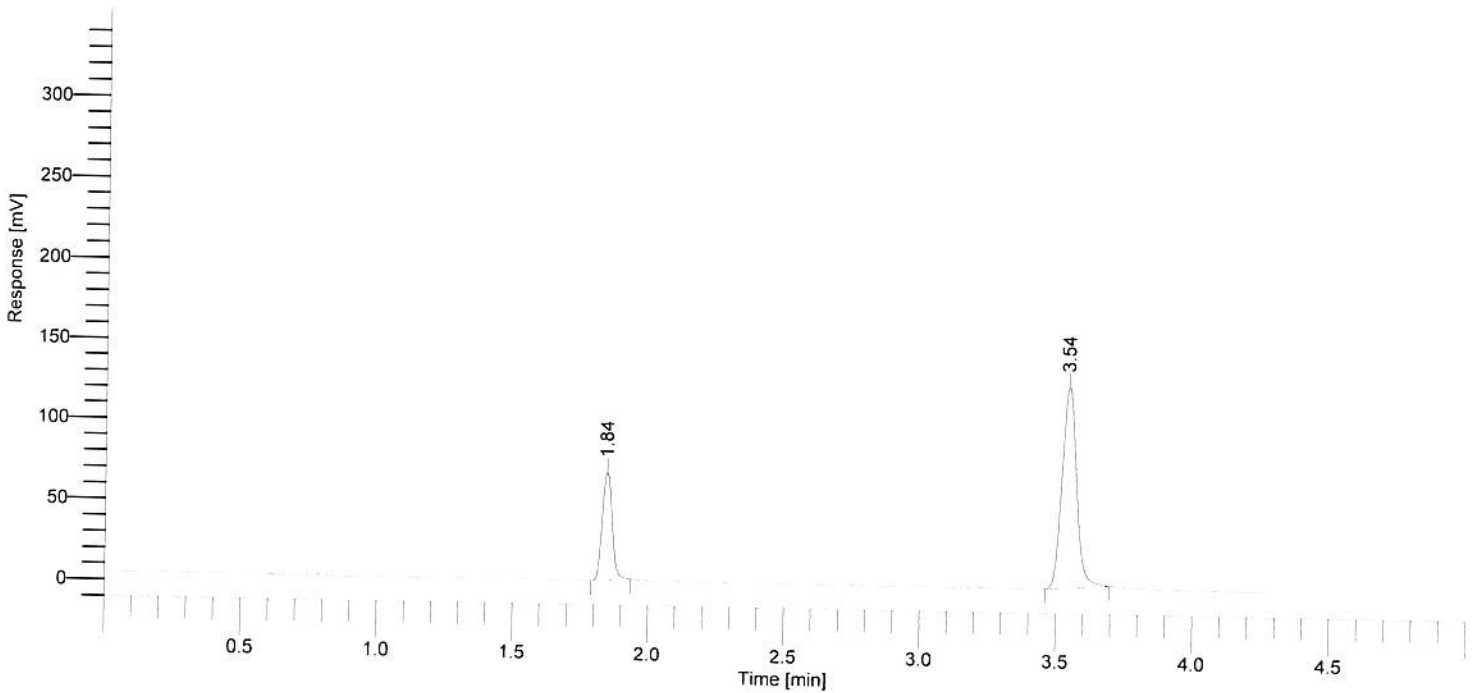

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.84
N-PROPANOL	3.54

15 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33286
Sample Name : GUTH 17310 #6
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 7

Date : 10/5/2017 3:13:14 PM
Data Acquisition Time : 10/5/2017 3:08:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-007.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

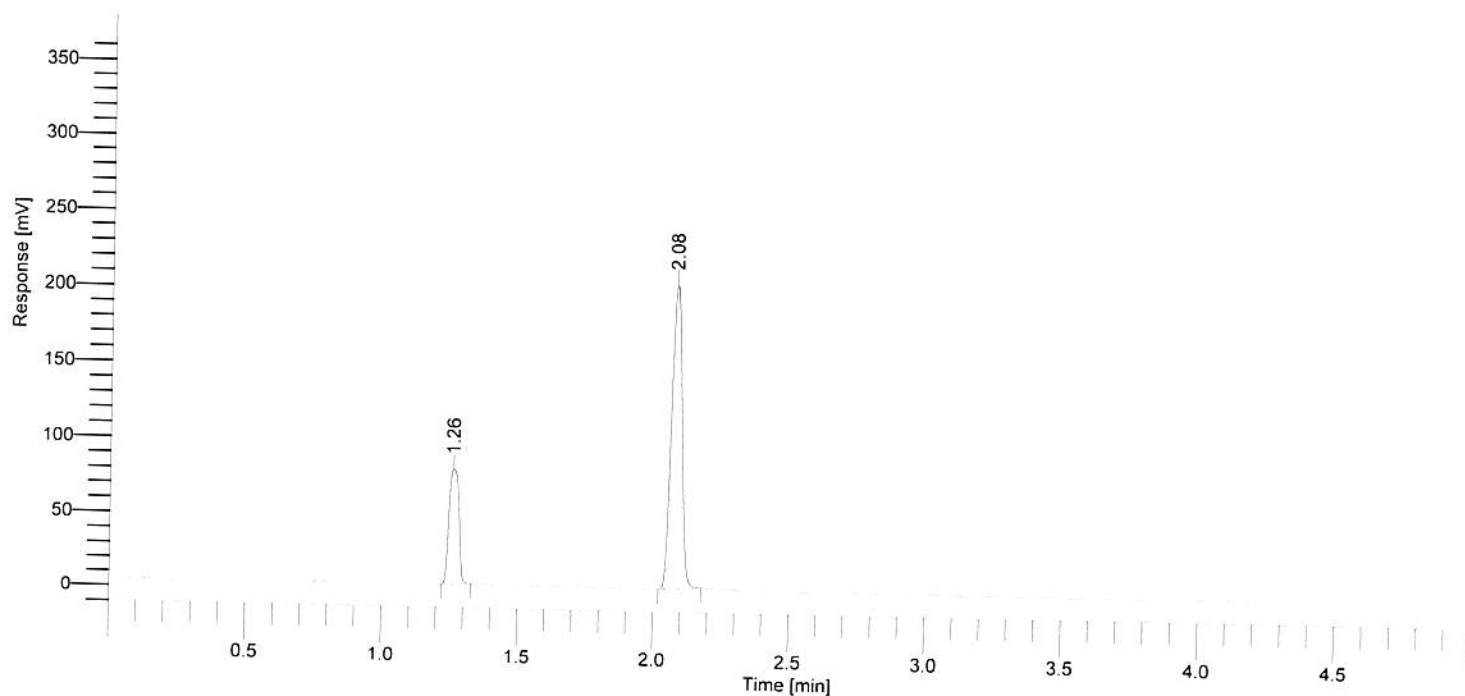

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1188	1.26	192774.0
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	541707.7

16 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33287
Sample Name : GUTH 17310 #6
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 7

Date : 10/5/2017 3:13:15 PM
Data Acquisition Time : 10/5/2017 3:08:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-007.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

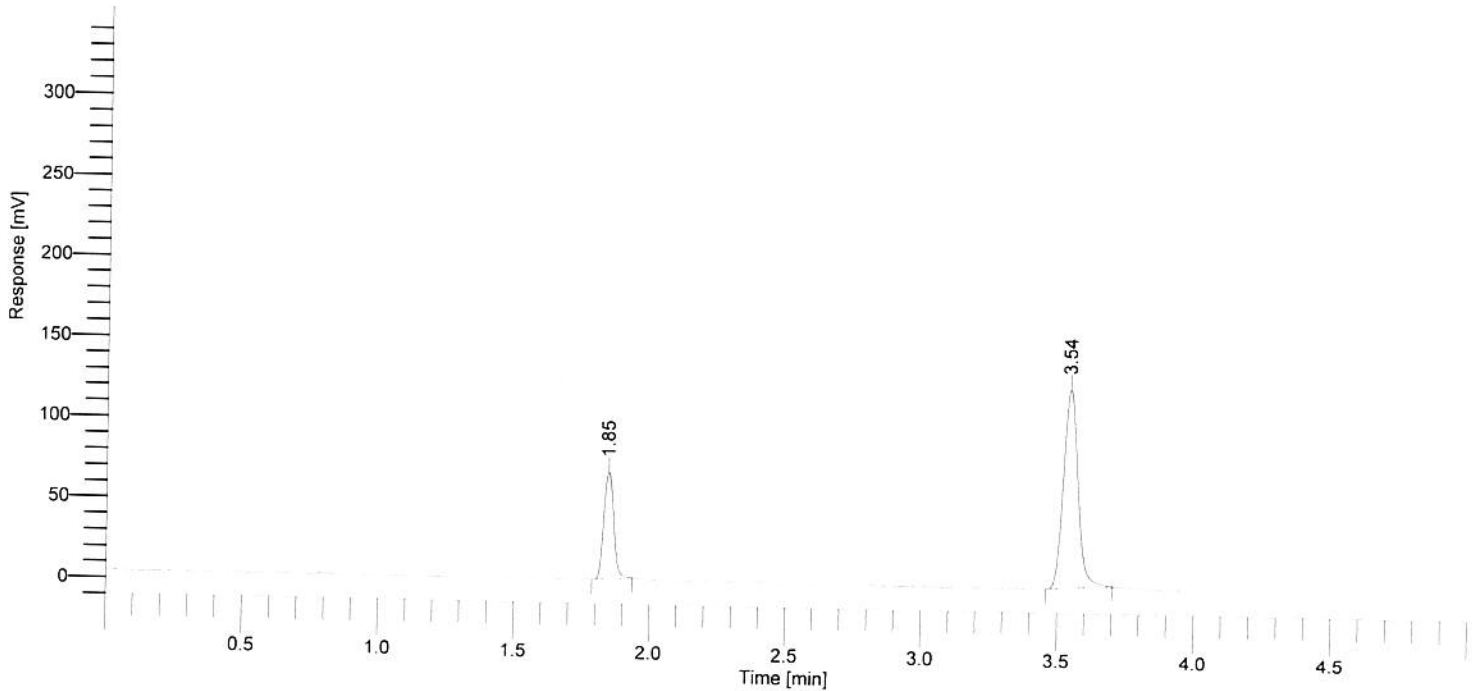

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

17 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33288
Sample Name : GUTH 17310 #7
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 8

Date : 10/5/2017 3:18:14 PM
Data Acquisition Time : 10/5/2017 3:13:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-008.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

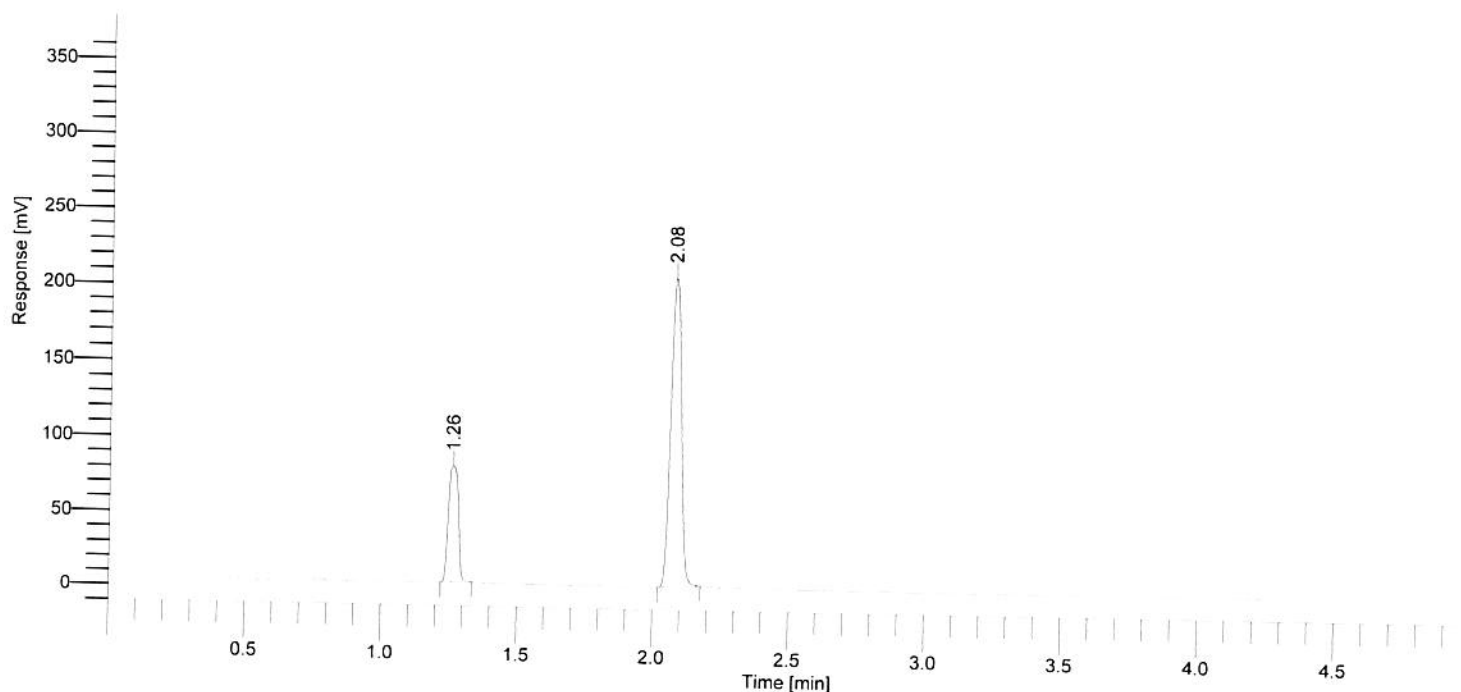

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1188	1.26	196115.6
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	551173.5

18 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33289
Sample Name : GUTH 17310 #7
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 8

Date : 10/5/2017 3:18:15 PM
Data Acquisition Time : 10/5/2017 3:13:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-008.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

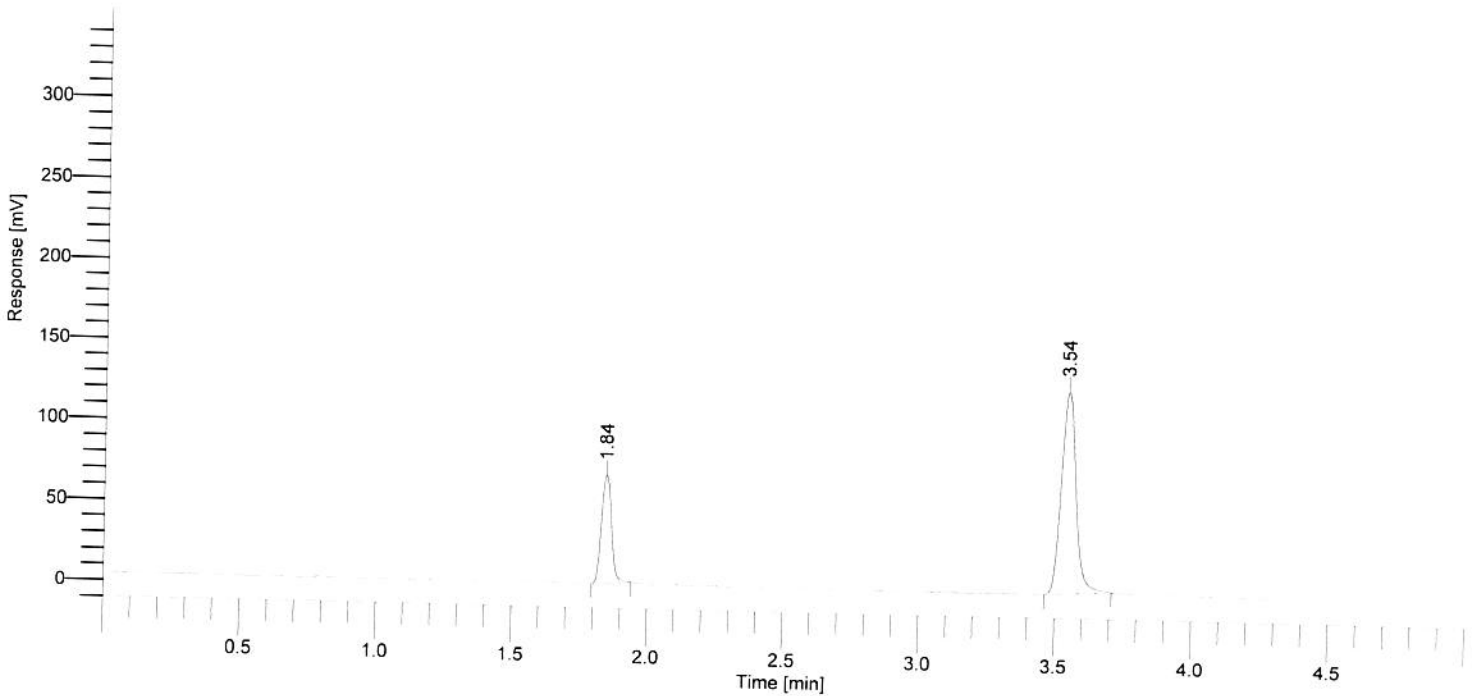

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.84
N-PROPANOL	3.54

19C

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33290
Sample Name : GUTH 17310 #7
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 9

Date : 10/5/2017 3:23:13 PM
Data Acquisition Time : 10/5/2017 3:18:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-009.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

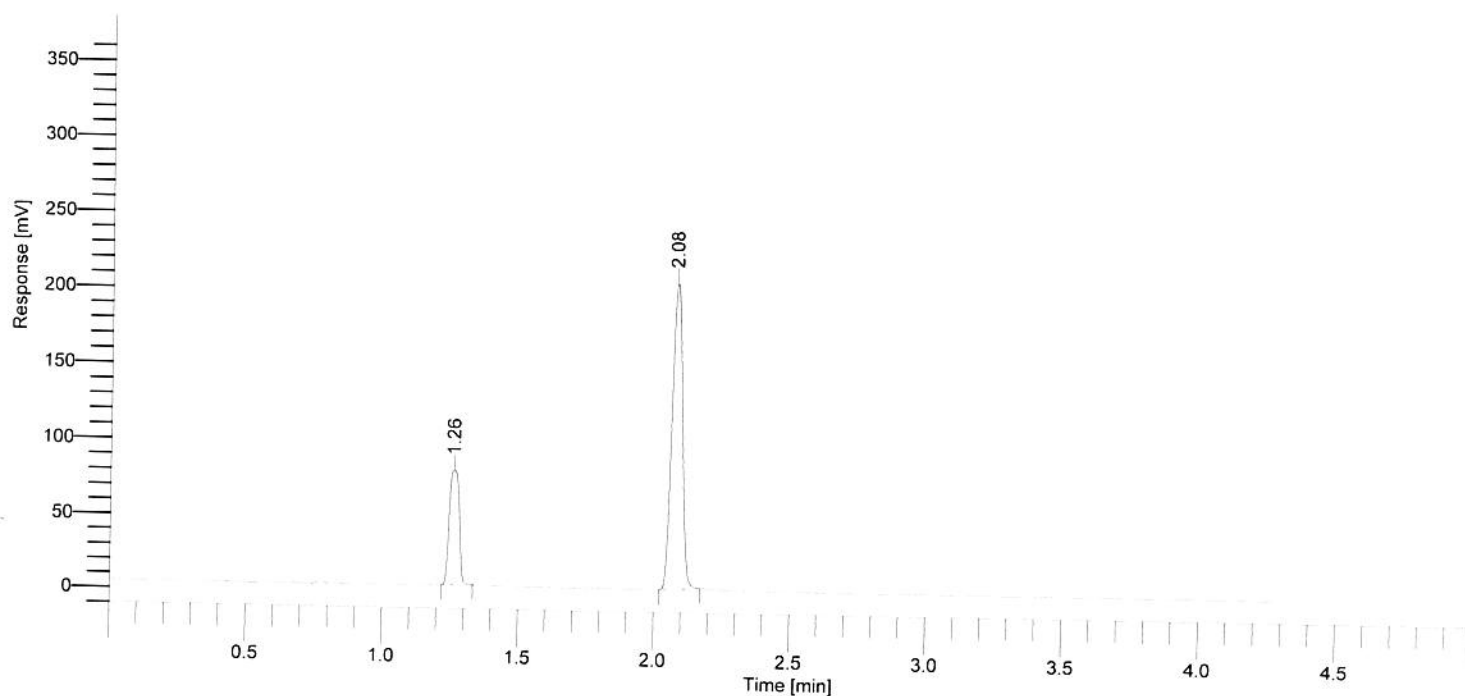

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1188	1.26	193915.8
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	544655.0

20 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33291
Sample Name : GUTH 17310 #7
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 9

Date : 10/5/2017 3:23:15 PM
Data Acquisition Time : 10/5/2017 3:18:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-009.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

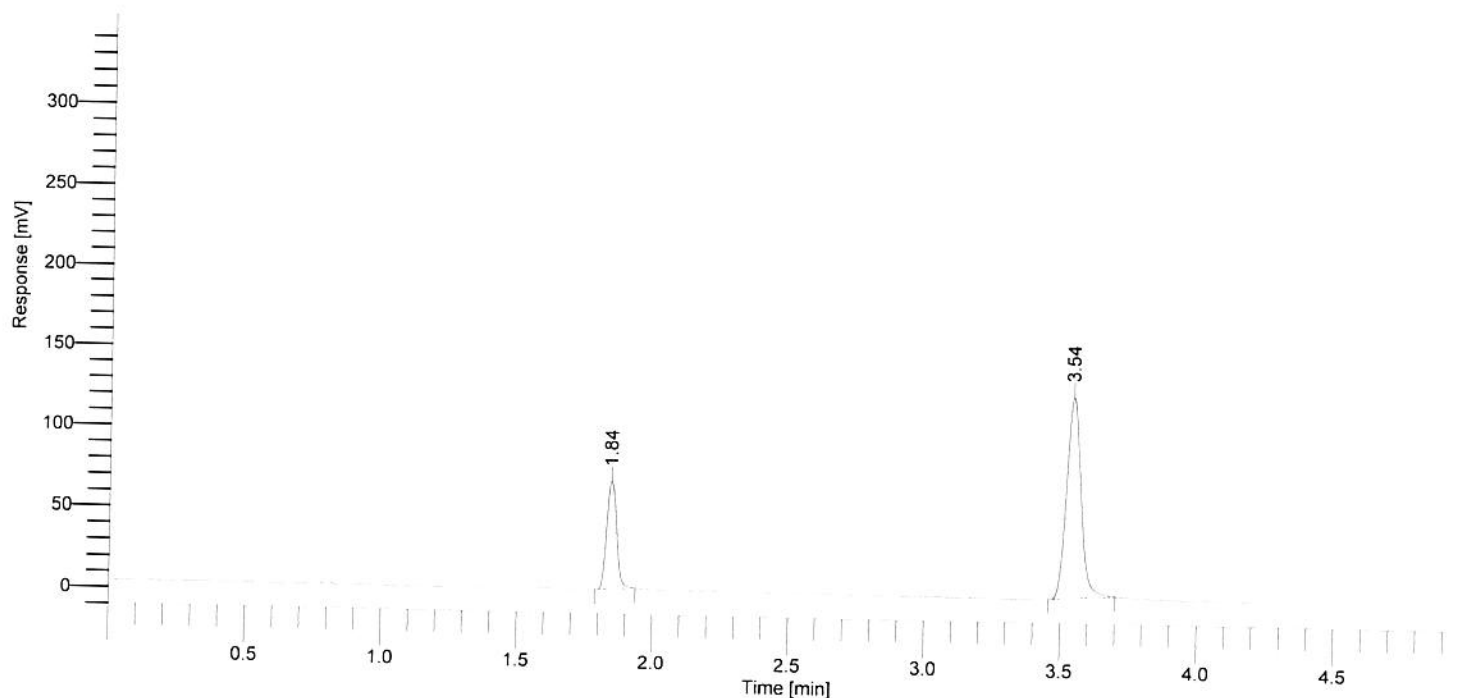

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.84
N-PROPANOL	3.54

21 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33292
Sample Name : GUTH 17310 #7
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 10

Date : 10/5/2017 3:28:14 PM
Data Acquisition Time : 10/5/2017 3:23:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-010.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

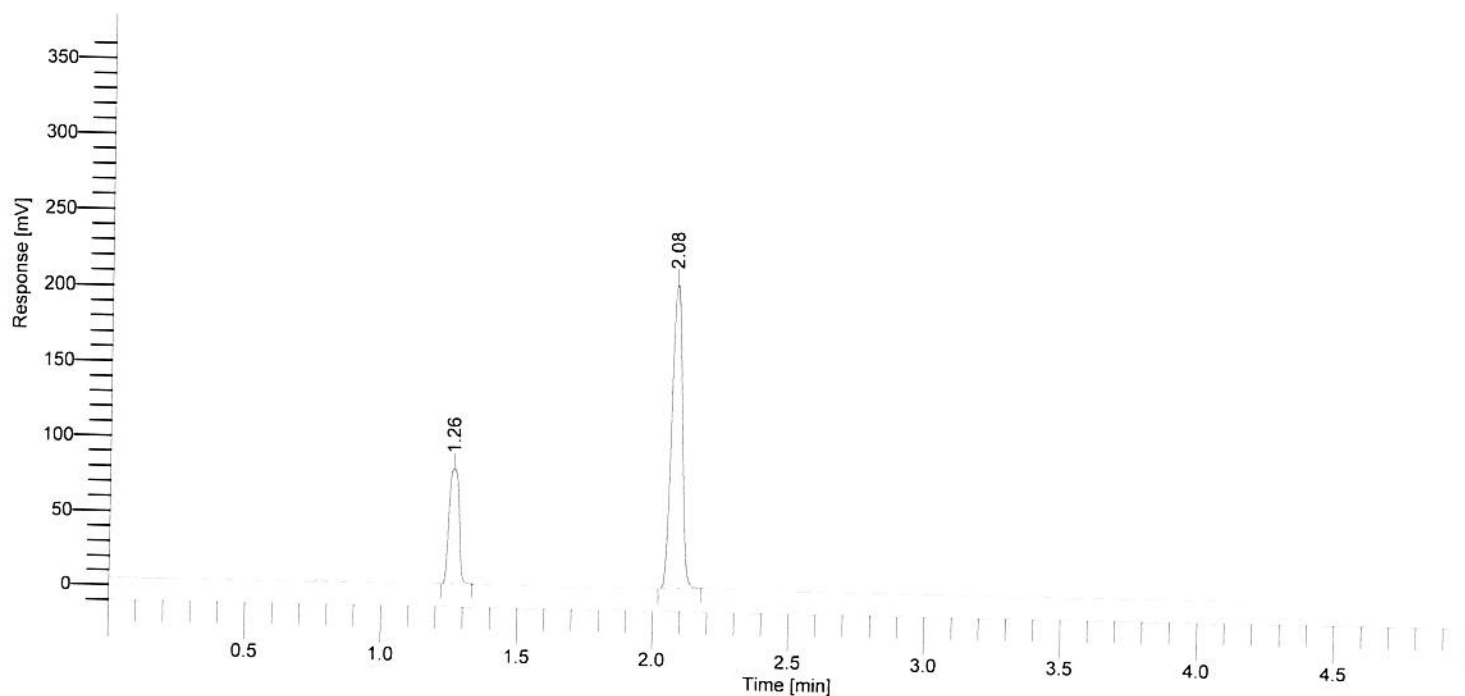

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1187	1.26	194348.1
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	546648.6

22 *CS*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33293
Sample Name : GUTH 17310 #7
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 10

Date : 10/5/2017 3:28:15 PM
Data Acquisition Time : 10/5/2017 3:23:48 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-010.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

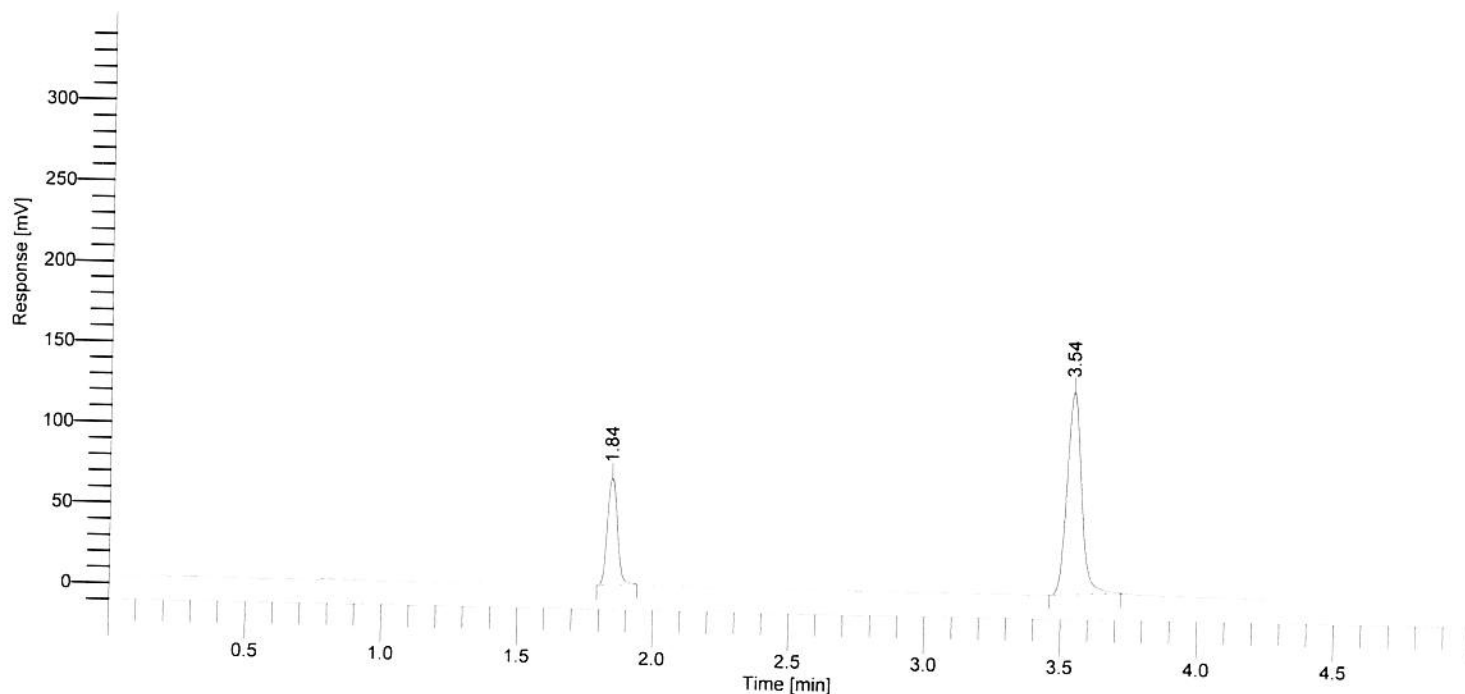

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.84
N-PROPANOL	3.54

23

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33294
Sample Name : GUTH 17310 #705
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 11

Date : 10/5/2017 3:33:13 PM
Data Acquisition Time : 10/5/2017 3:28:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-011.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

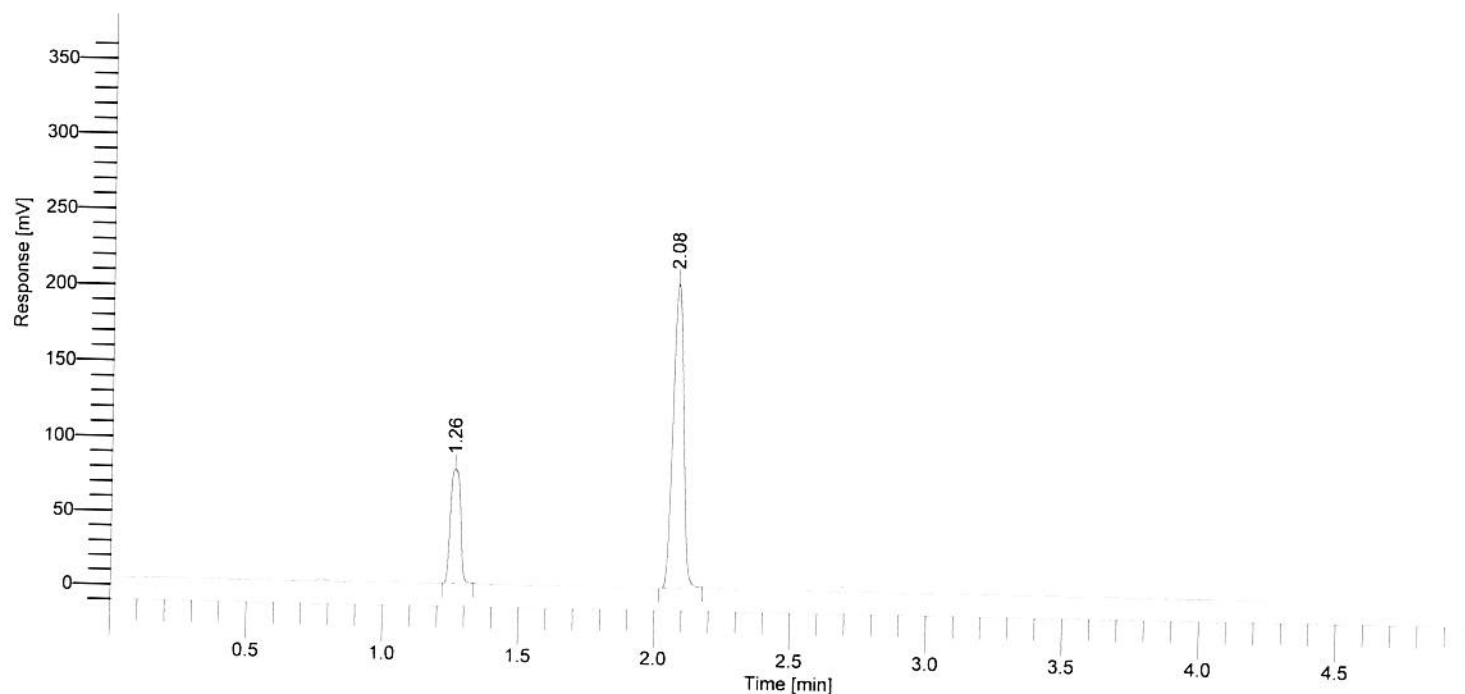

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1193	1.26	193164.7
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	540597.8

24 cm

Software Version : 6.3.2.0646
 Reprocess Number : d-2ua3120y7b: 33295
 Sample Name : GUTH 17310 #705
 Instrument Name : GCHeadspace2
 Rack/Vial : 0/0
 Sample Amount : 1.000000
 Cycle : 11

Date : 10/5/2017 3:33:15 PM
 Data Acquisition Time : 10/5/2017 3:28:48 PM
 Channel : B
 Operator : manager
 Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-011.rst
 Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

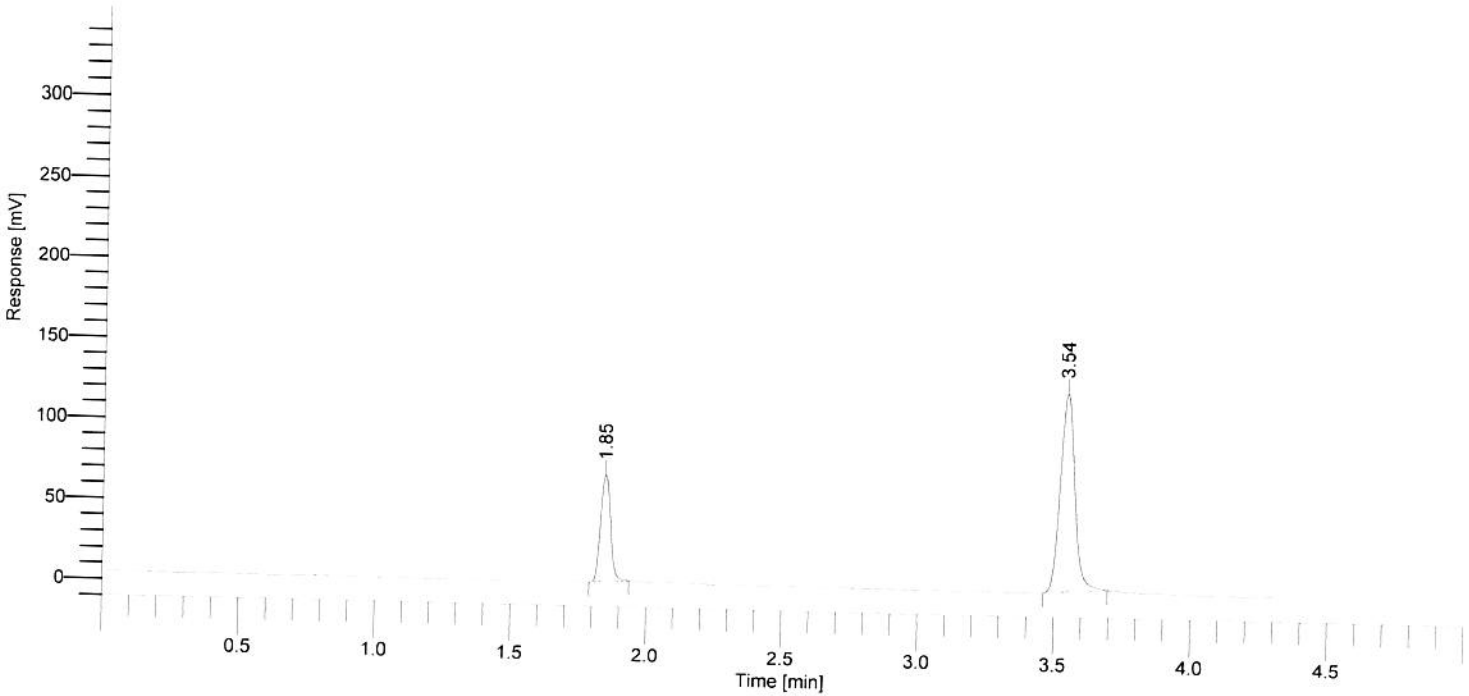

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

25 Gm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33296
Sample Name : GUTH 17310 #705
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 12

Date : 10/5/2017 3:38:14 PM
Data Acquisition Time : 10/5/2017 3:33:48 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-012.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

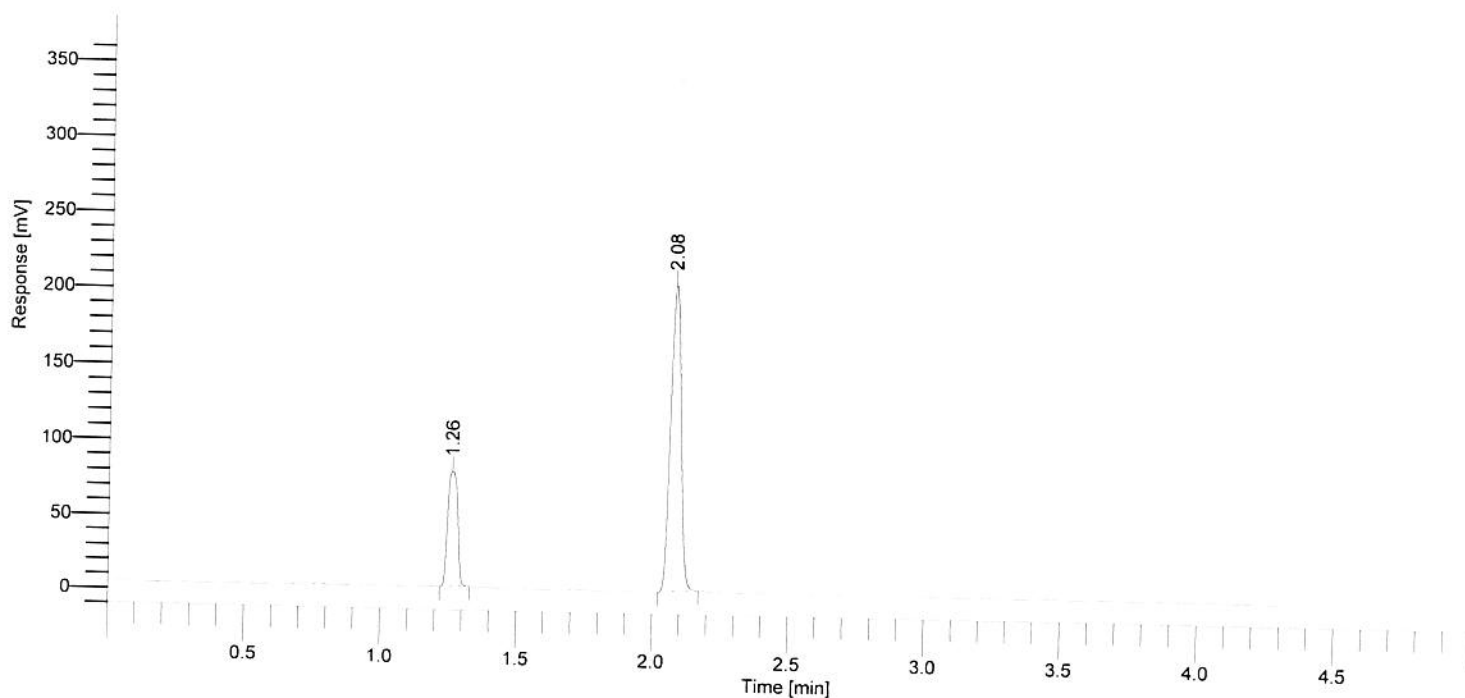

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1190	1.26	193190.0
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	541675.7

266

Software Version : 6.3.2.0646
 Reprocess Number : d-2ua3120y7b: 33297
 Sample Name : GUTH 17310 #705
 Instrument Name : GCHeadspace2
 Rack/Vial : 0/0
 Sample Amount : 1.000000
 Cycle : 12

Date : 10/5/2017 3:38:15 PM
 Data Acquisition Time : 10/5/2017 3:33:48 PM
 Channel : B
 Operator : manager
 Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-012.rst
 Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

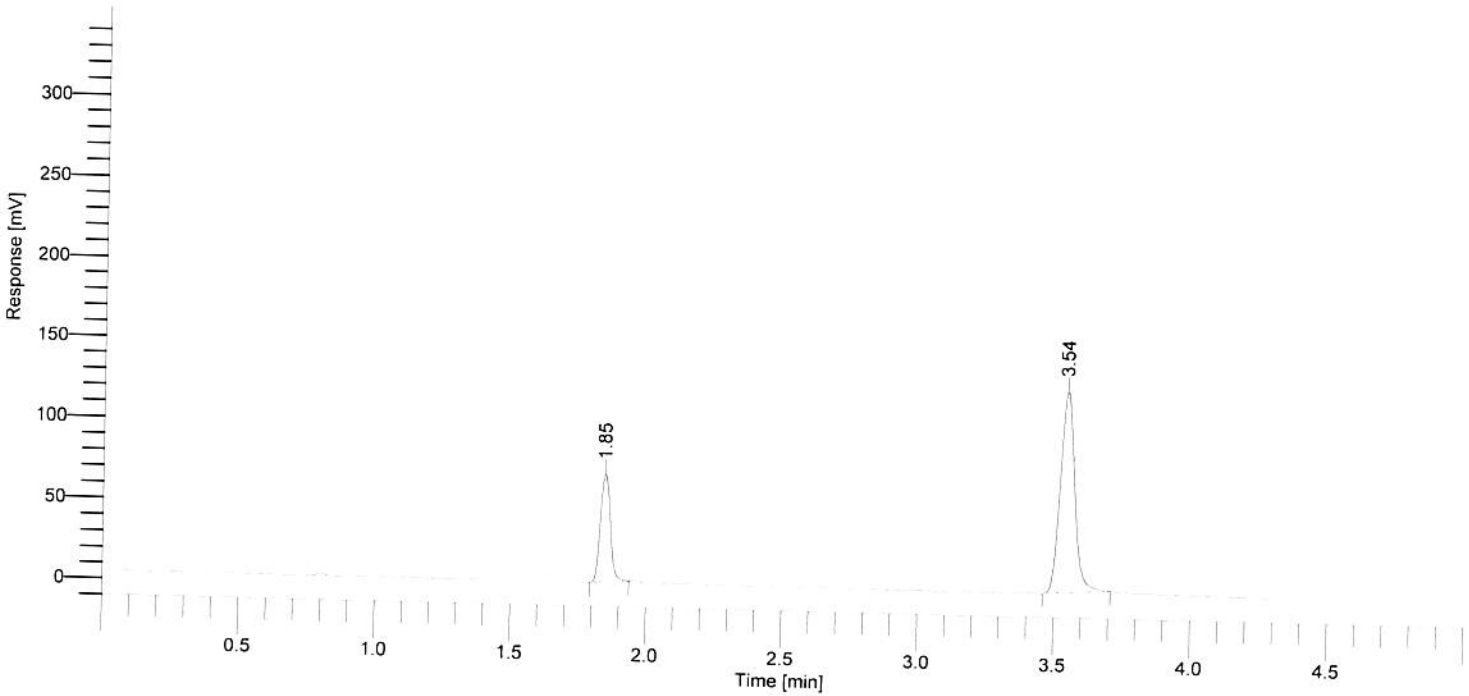

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

27 *[Signature]*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33298
Sample Name : GUTH 17310 #705
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 13

Date : 10/5/2017 3:43:13 PM
Data Acquisition Time : 10/5/2017 3:38:49 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-013.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

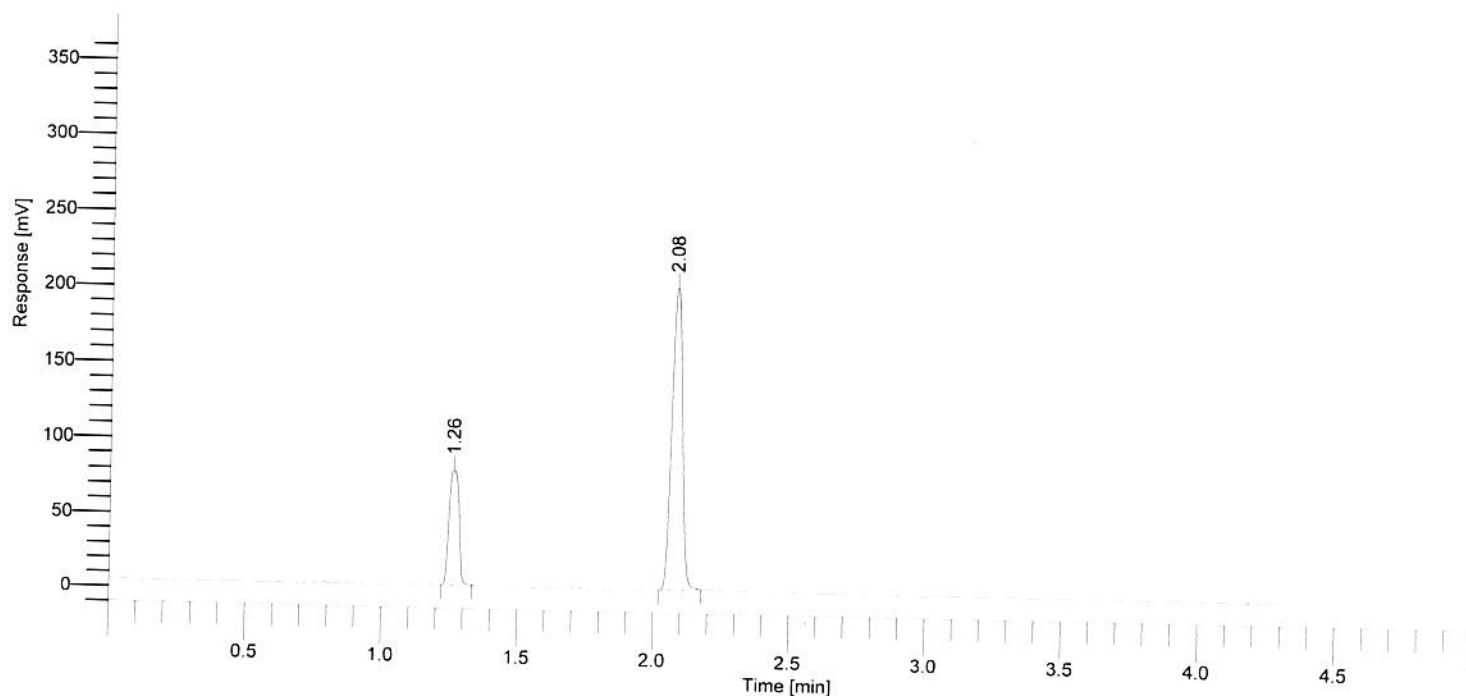

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1191	1.26	192690.0
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	539970.0

28 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33299
Sample Name : GUTH 17310 #705
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 13

Date : 10/5/2017 3:43:15 PM
Data Acquisition Time : 10/5/2017 3:38:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-013.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

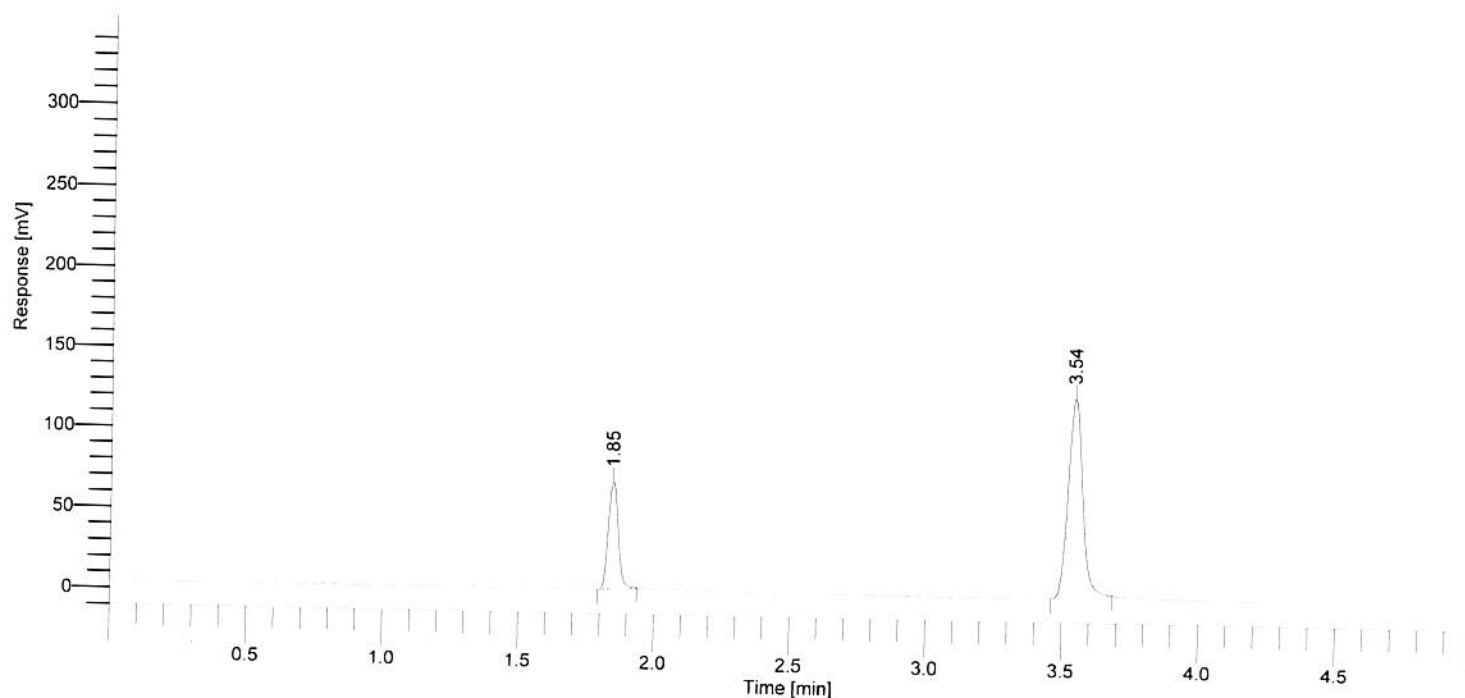

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

29 *[Signature]*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33300
Sample Name : 0.10% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 14

Date : 10/5/2017 3:48:14 PM
Data Acquisition Time : 10/5/2017 3:43:49 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-014.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

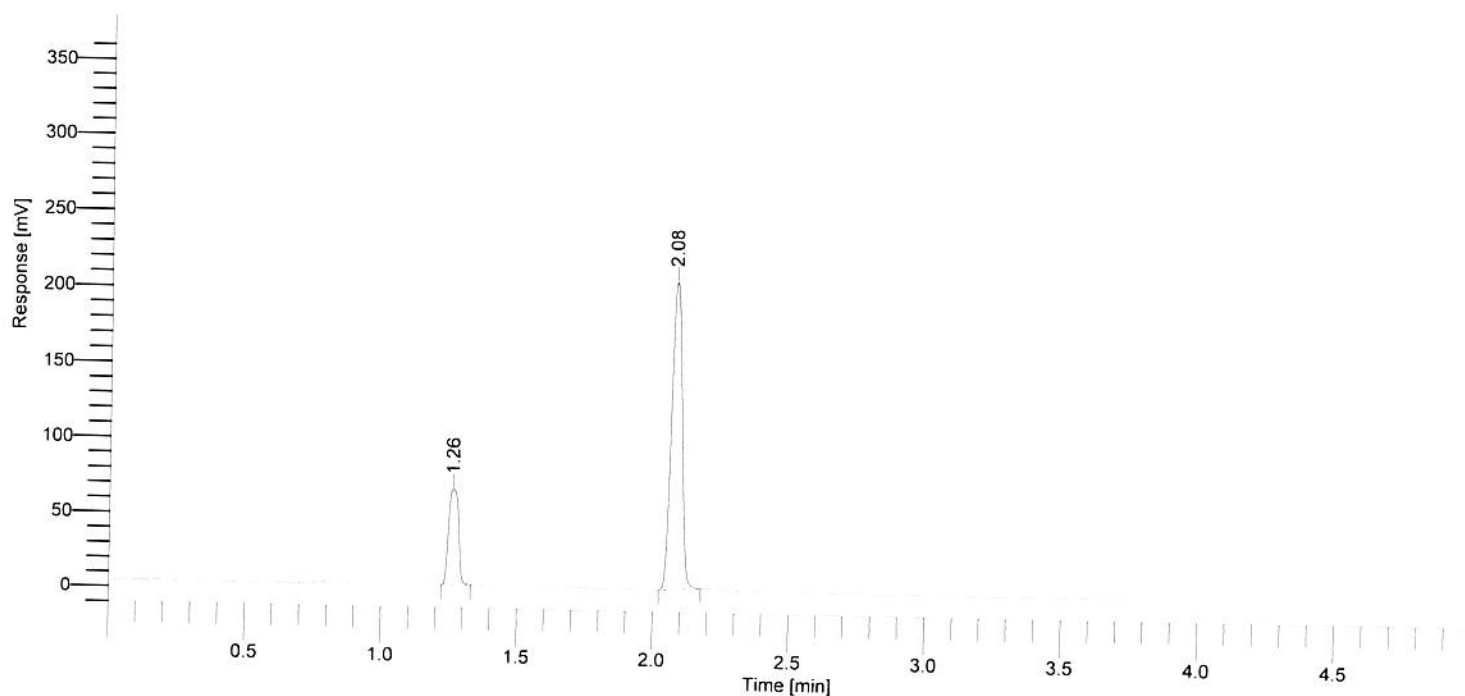

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.0978	1.26	160276.3
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	549078.5

30 *cm*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33301
Sample Name : 0.10% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 14

Date : 10/5/2017 3:48:15 PM
Data Acquisition Time : 10/5/2017 3:43:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-014.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

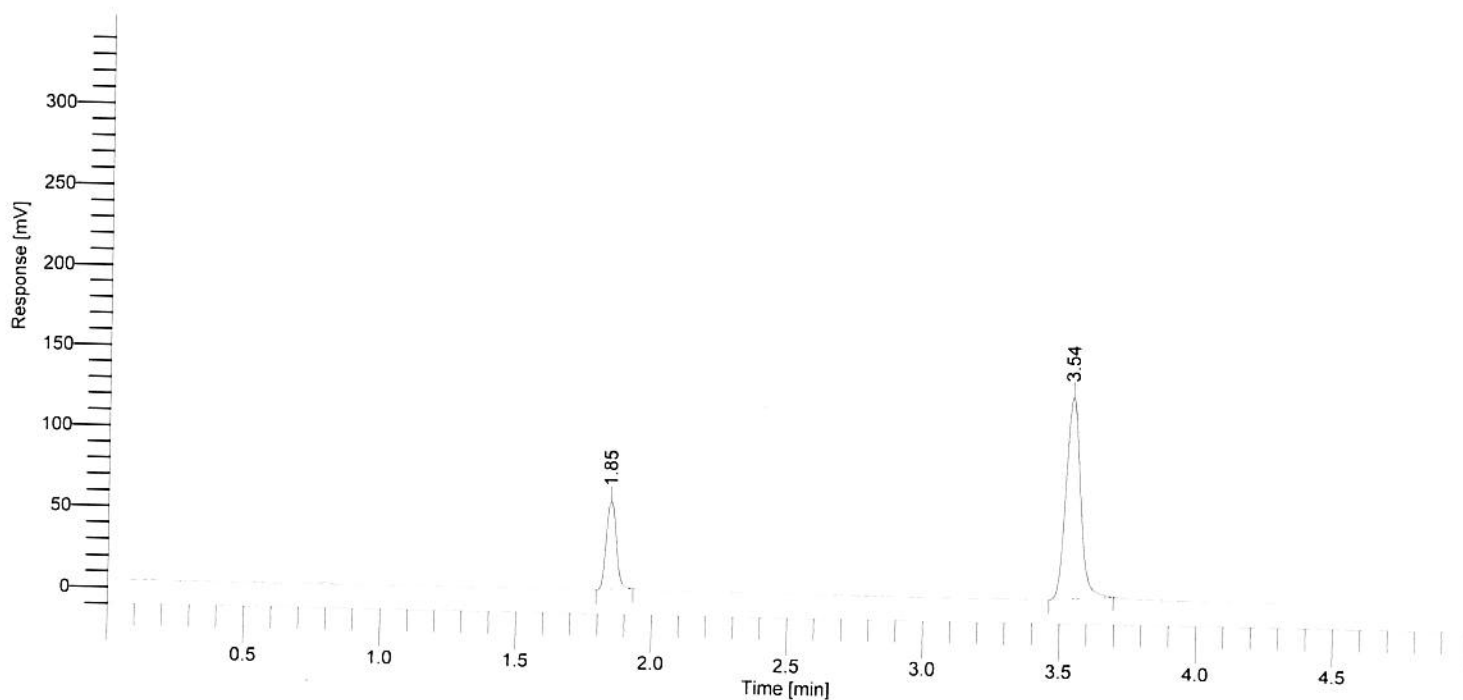

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

31

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33302
Sample Name : NEG CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 15

Date : 10/5/2017 3:53:13 PM
Data Acquisition Time : 10/5/2017 3:48:49 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-015.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

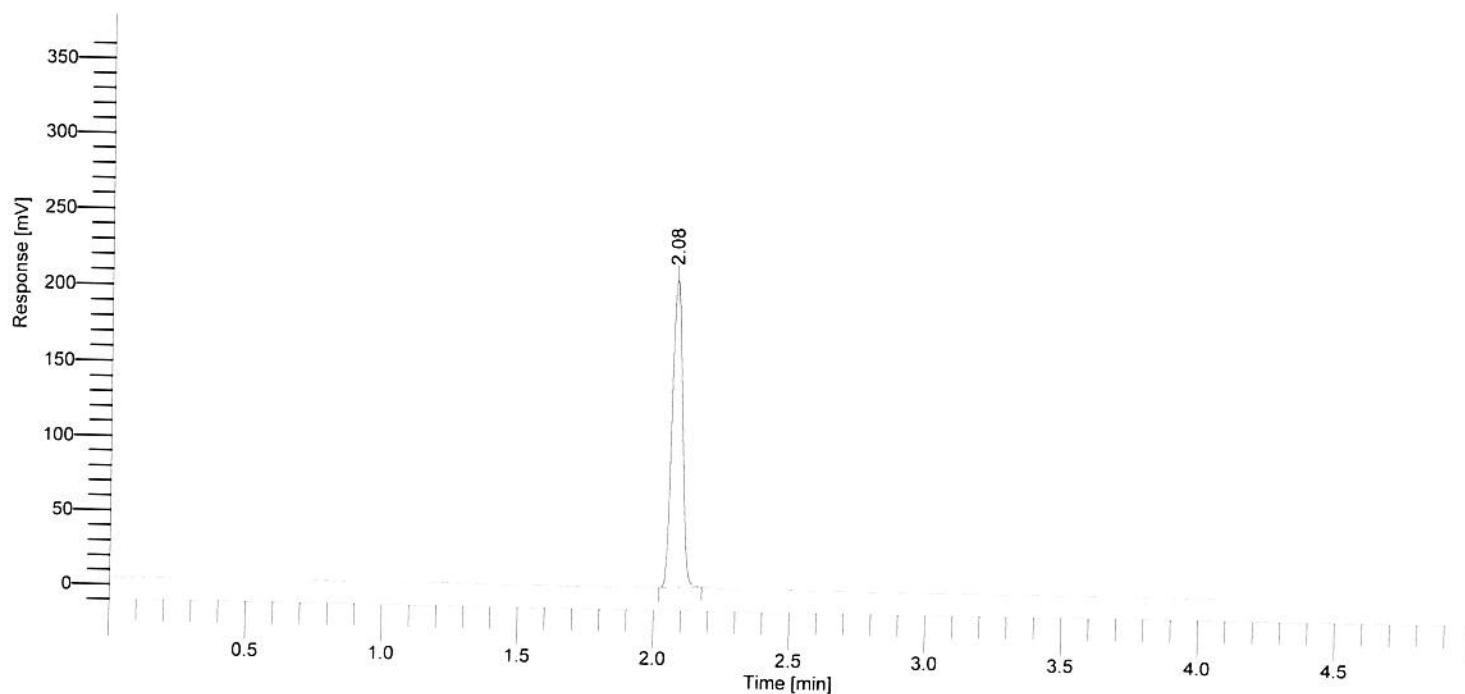

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.0000	1.26	0.0
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	546695.4

32 *CSM*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33303
Sample Name : NEG CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 15

Date : 10/5/2017 3:53:15 PM
Data Acquisition Time : 10/5/2017 3:48:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-015.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

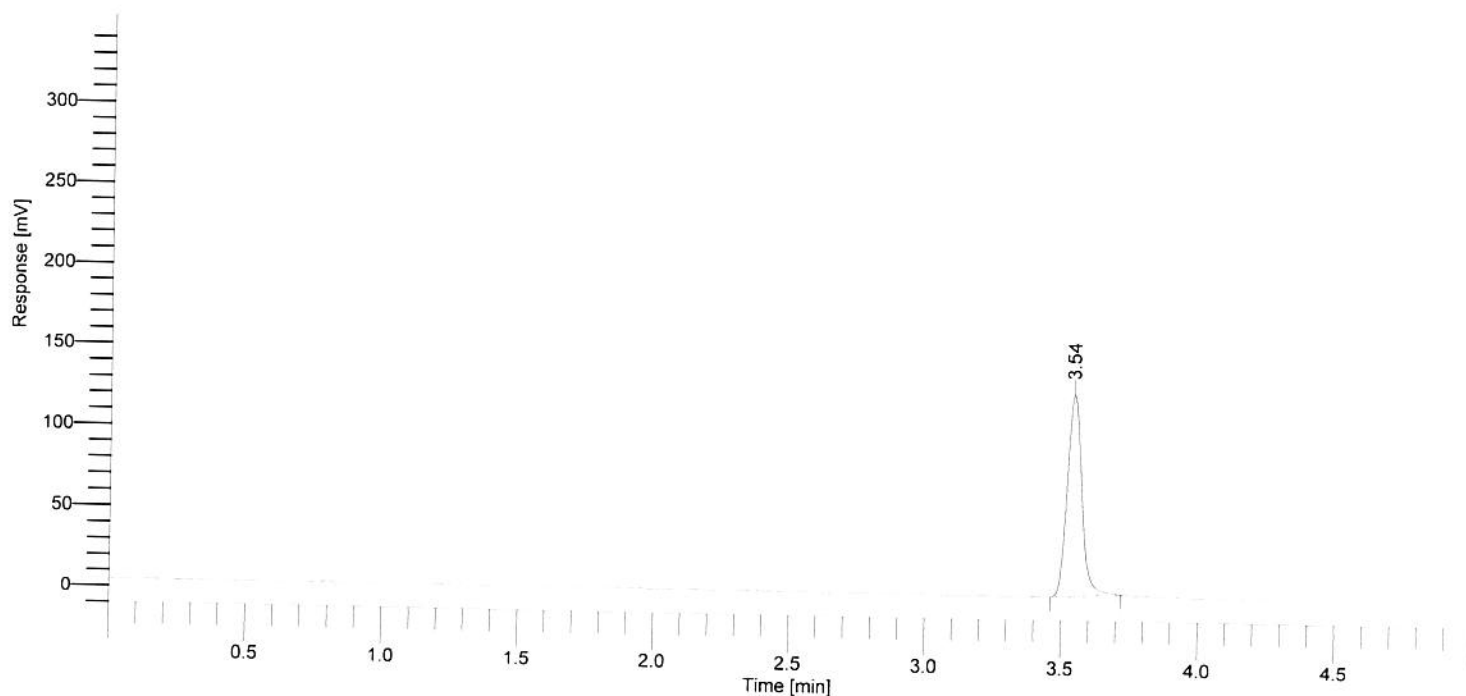

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
N-PROPANOL	3.54

33 *CSM*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33304
Sample Name : GUTH 17310 #706
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 16

Date : 10/5/2017 3:58:14 PM
Data Acquisition Time : 10/5/2017 3:53:49 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-016.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

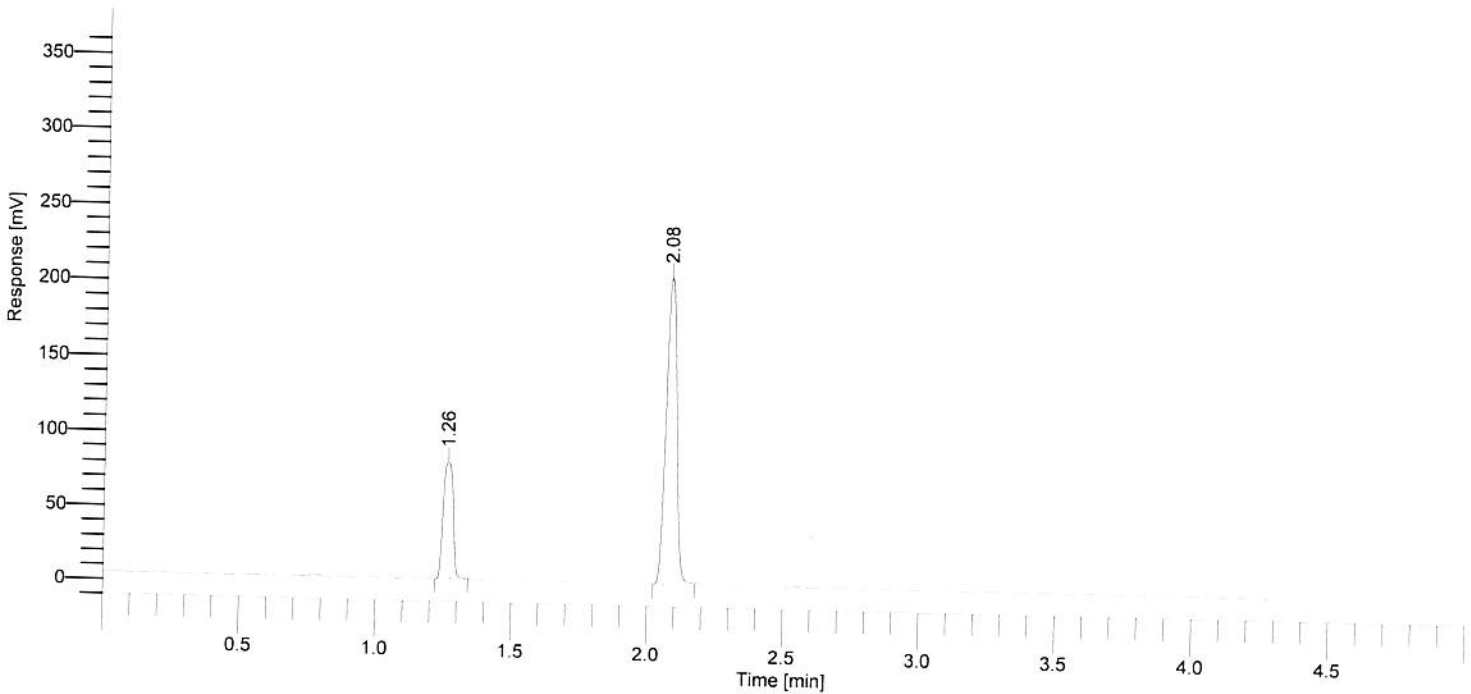

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1191	1.26	195761.5
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	548750.8

34 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33305
Sample Name : GUTH 17310 #706
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 16

Date : 10/5/2017 3:58:15 PM
Data Acquisition Time : 10/5/2017 3:53:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-016.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

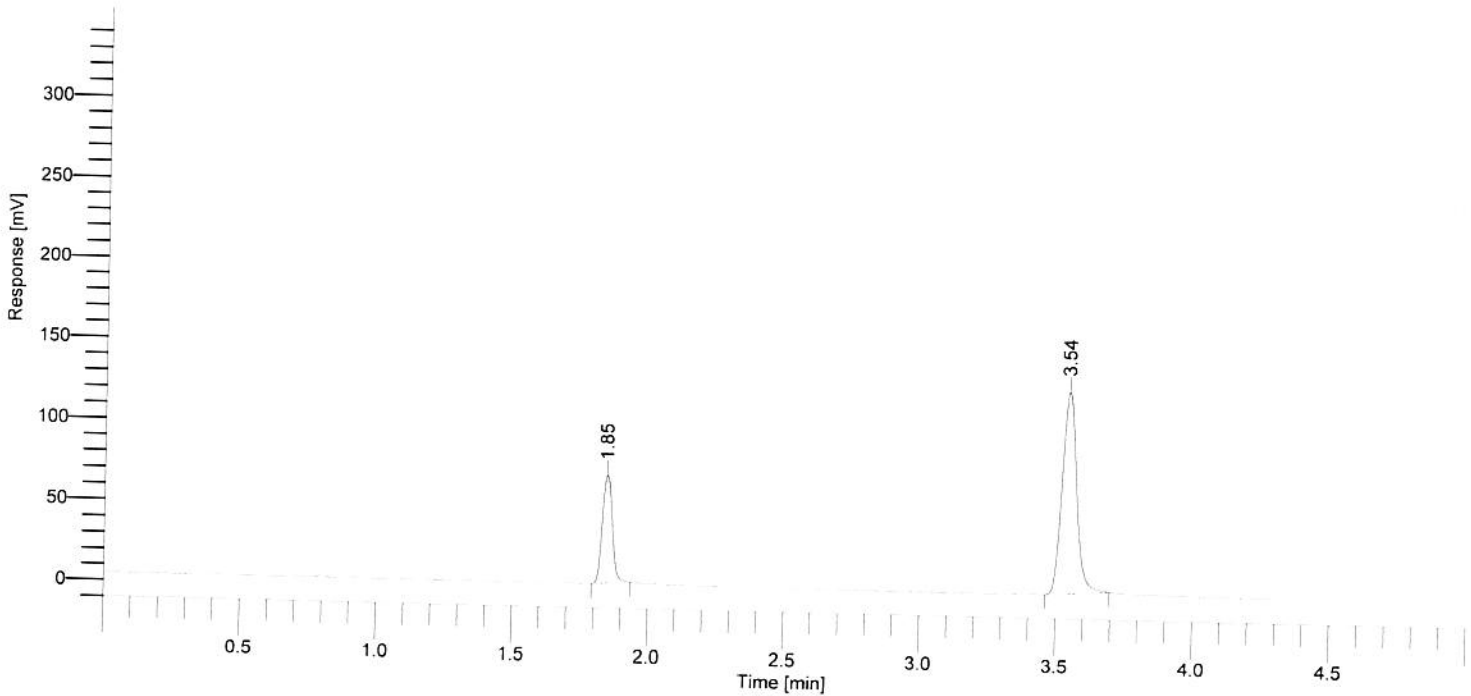

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

35 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33306
Sample Name : GUTH 17310 #706
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 17

Date : 10/5/2017 4:03:13 PM
Data Acquisition Time : 10/5/2017 3:58:49 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-017.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

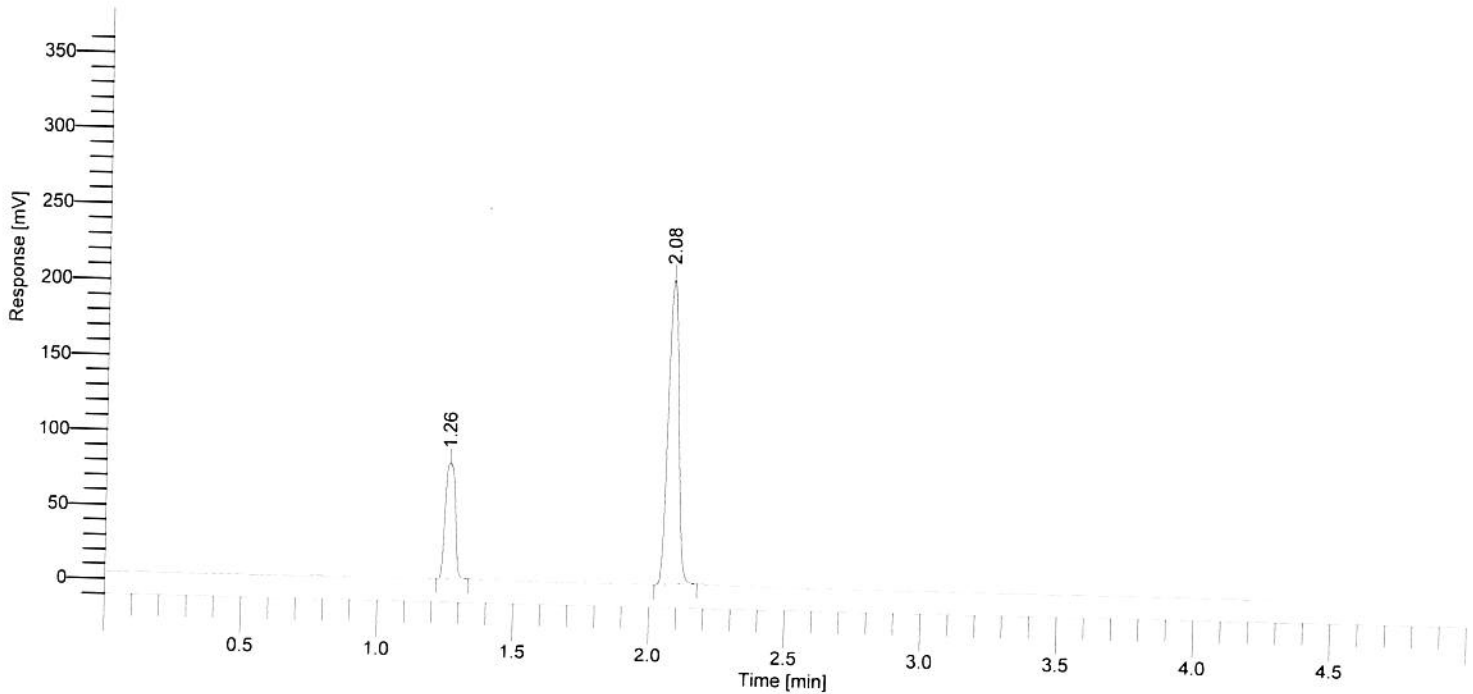

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1190	1.26	194109.5
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	544277.4

36 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33307
Sample Name : GUTH 17310 #706
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 17

Date : 10/5/2017 4:03:15 PM
Data Acquisition Time : 10/5/2017 3:58:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-017.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

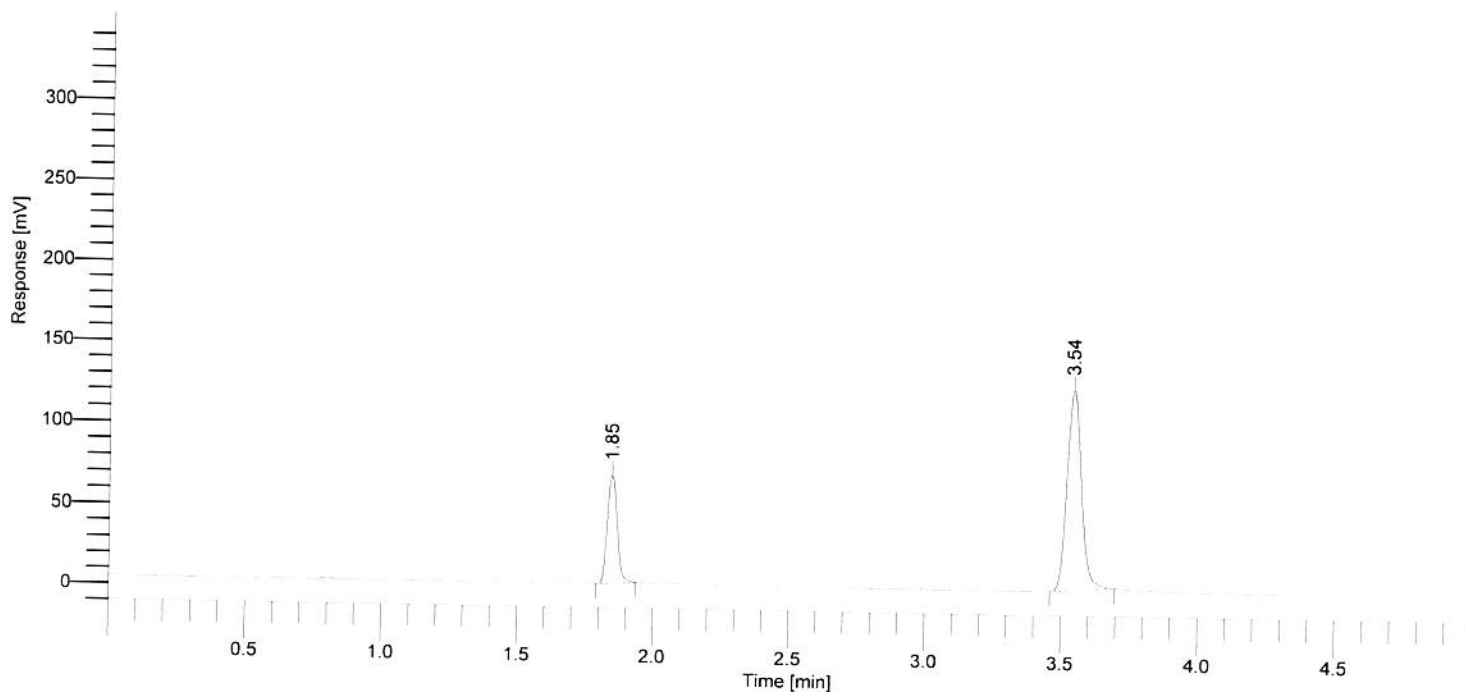

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

37 *CSM*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33308
Sample Name : GUTH 17310 #706
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 18

Date : 10/5/2017 4:08:13 PM
Data Acquisition Time : 10/5/2017 4:03:49 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-018.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

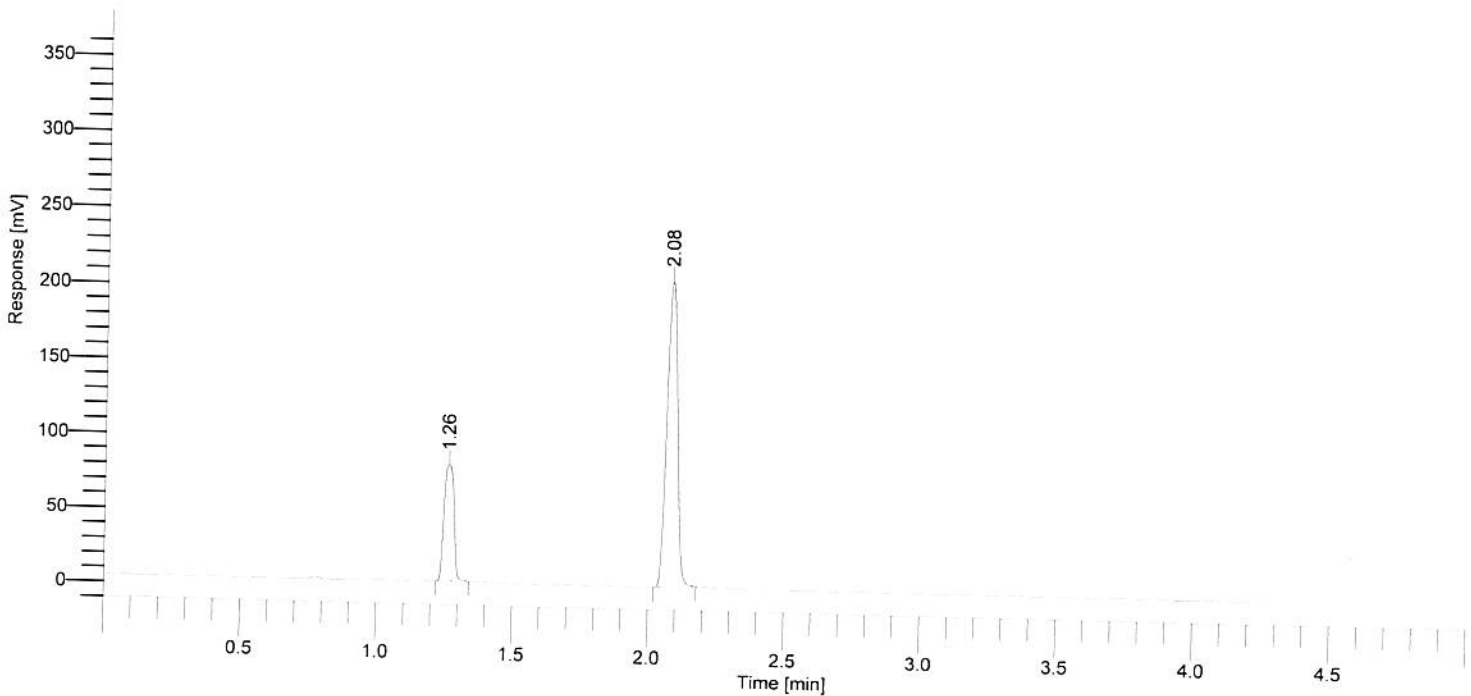

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1191	1.26	195093.9
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	546878.2

38 *Cam*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33309
Sample Name : GUTH 17310 #706
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 18

Date : 10/5/2017 4:08:15 PM
Data Acquisition Time : 10/5/2017 4:03:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-018.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

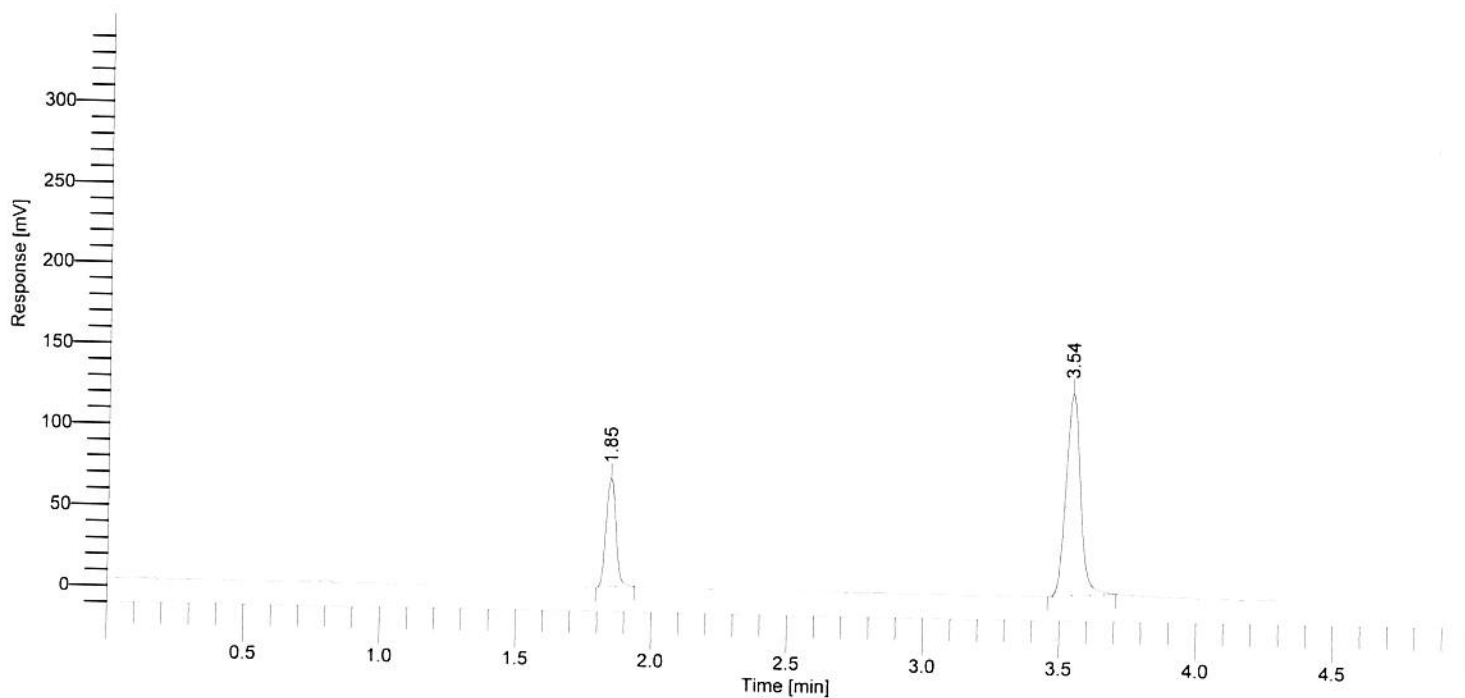

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

3912m

Software Version : 6.3.2.0646
 Reprocess Number : d-2ua3120y7b: 33310
 Sample Name : GUTH 17310 #1390
 Instrument Name : GCHeadspace2
 Rack/Vial : 0/0
 Sample Amount : 1.000000
 Cycle : 19

Date : 10/5/2017 4:13:13 PM
 Data Acquisition Time : 10/5/2017 4:08:49 PM
 Channel : A
 Operator : manager
 Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-019.rst
 Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

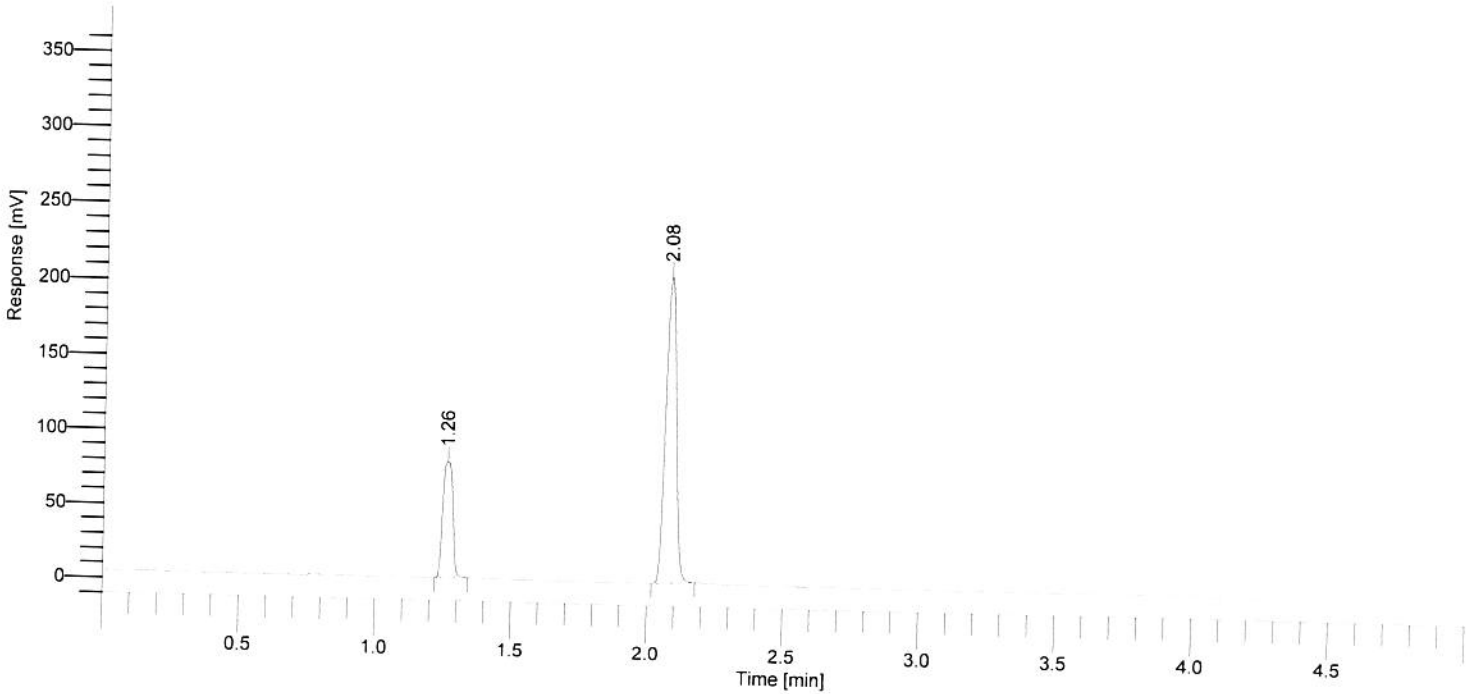

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μV·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1191	1.26	194872.1
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	545922.8

40 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33311
Sample Name : GUTH 17310 #1390
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 19

Date : 10/5/2017 4:13:15 PM
Data Acquisition Time : 10/5/2017 4:08:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-019.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

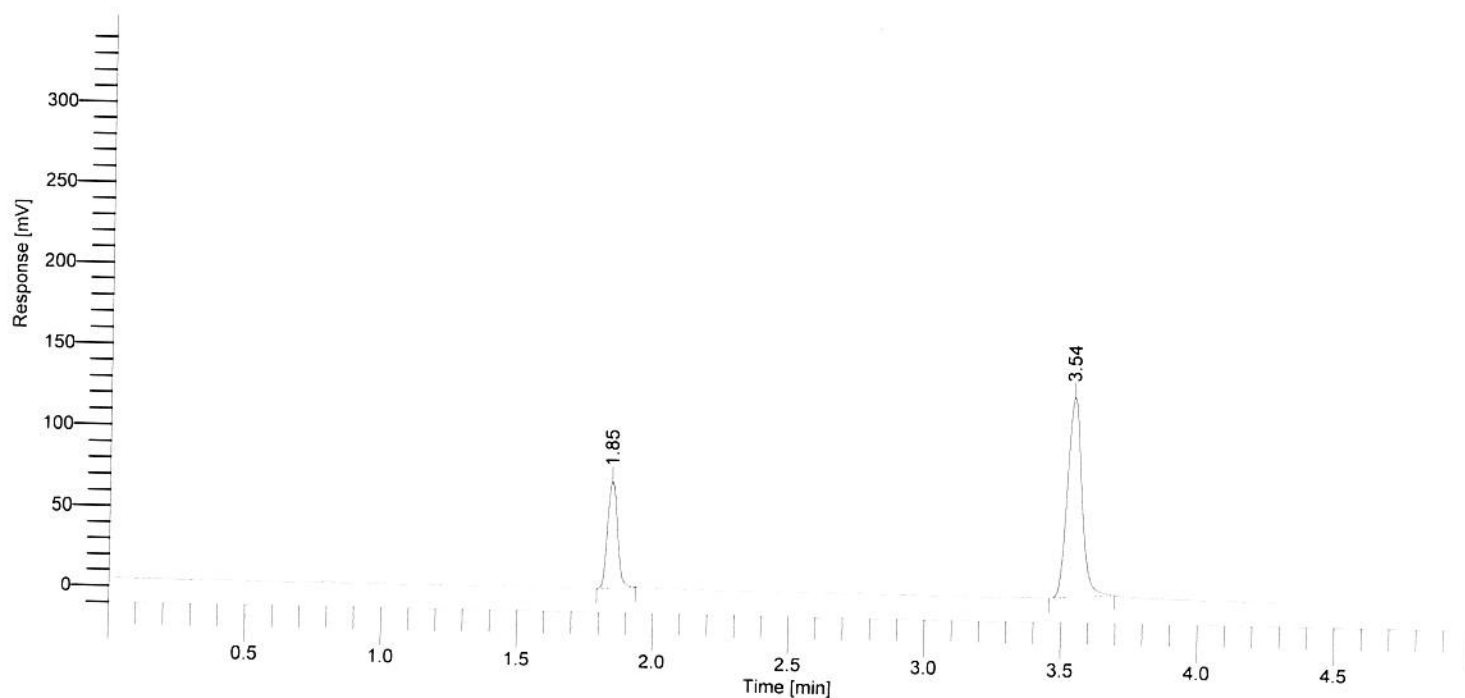

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

41Cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33312
Sample Name : GUTH 17310 #1390
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 20

Date : 10/5/2017 4:18:13 PM
Data Acquisition Time : 10/5/2017 4:13:49 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-020.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

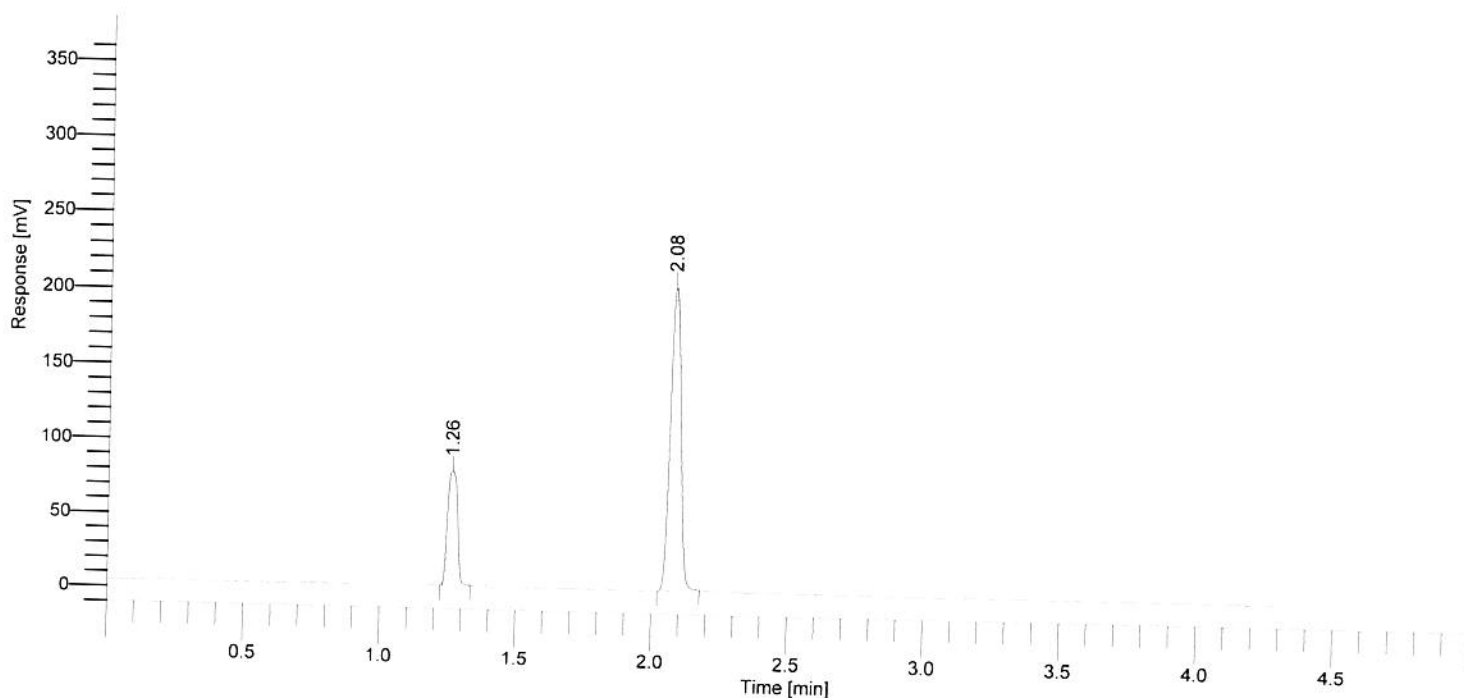

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1190	1.26	193212.5
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	542004.6

42 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33313
Sample Name : GUTH 17310 #1390
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 20

Date : 10/5/2017 4:18:15 PM
Data Acquisition Time : 10/5/2017 4:13:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-020.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

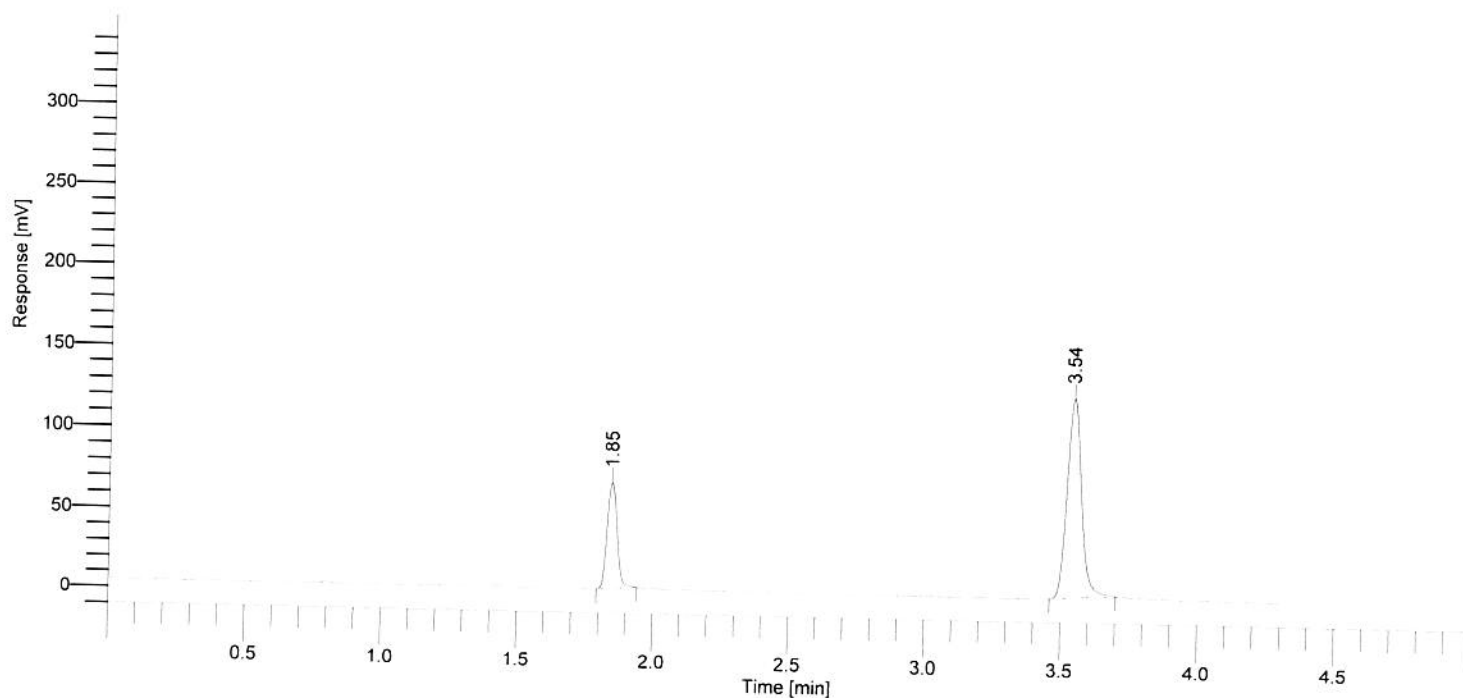

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

43 *Cam*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33314
Sample Name : GUTH 17310 #1390
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 21

Date : 10/5/2017 4:23:13 PM
Data Acquisition Time : 10/5/2017 4:18:49 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-021.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

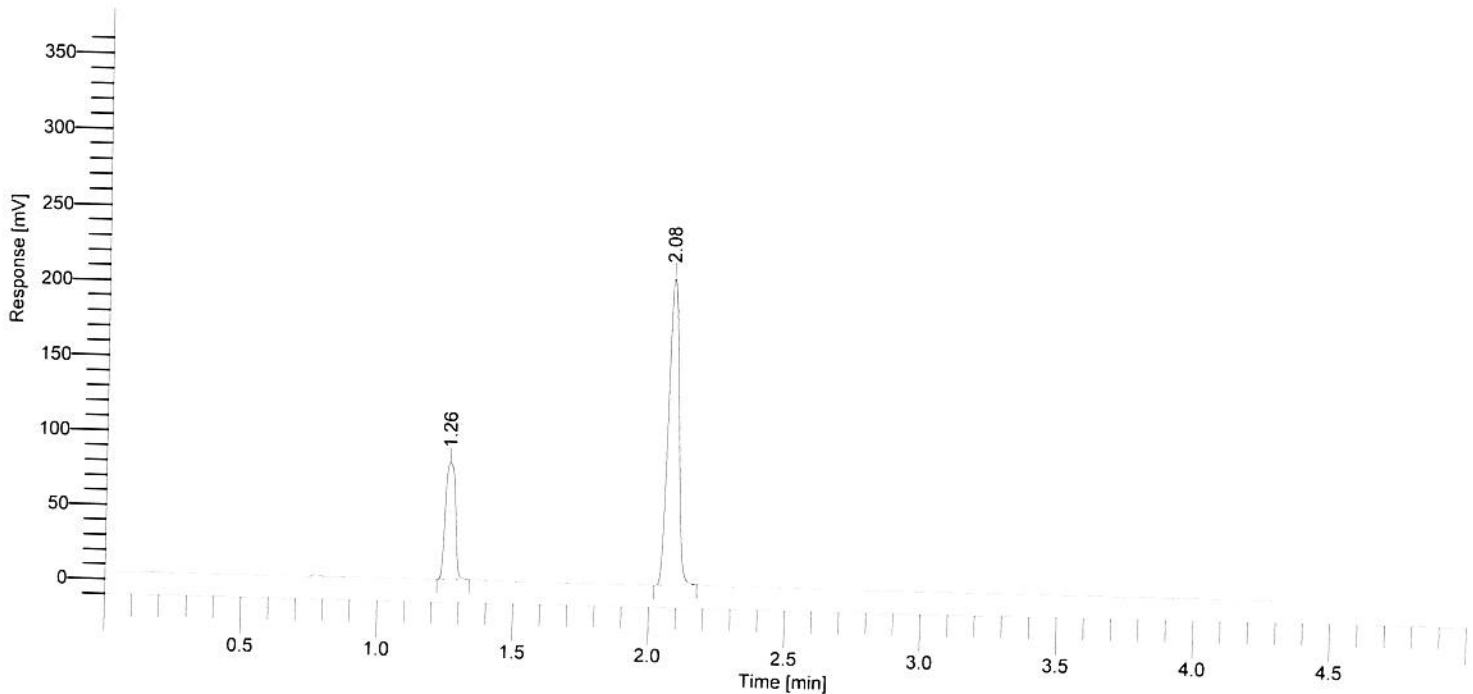

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1191	1.26	195906.7
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	548925.5

44 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33315
Sample Name : GUTH 17310 #1390
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 21

Date : 10/5/2017 4:23:15 PM
Data Acquisition Time : 10/5/2017 4:18:49 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-021.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

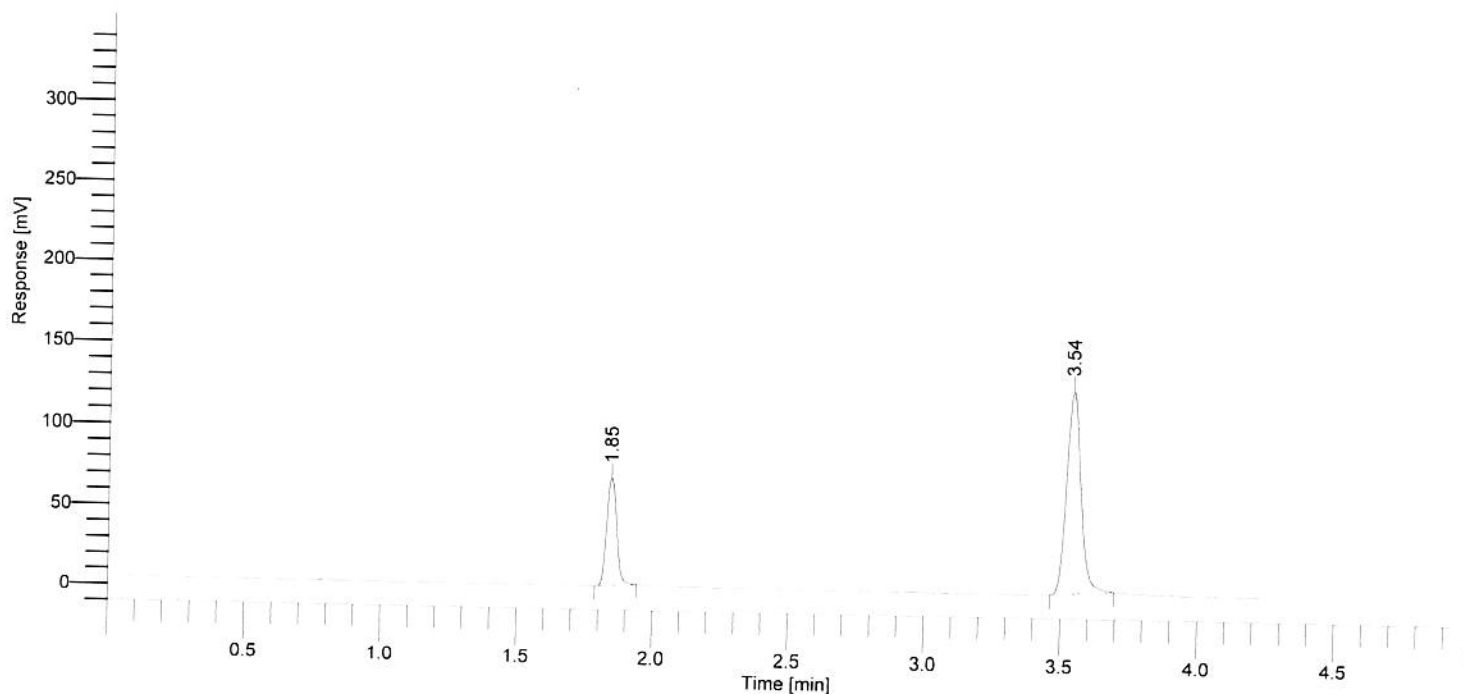

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

45 *CSM*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33316
Sample Name : GUTH 17310 #1391
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 22

Date : 10/5/2017 4:28:10 PM
Data Acquisition Time : 10/5/2017 4:23:46 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-022.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

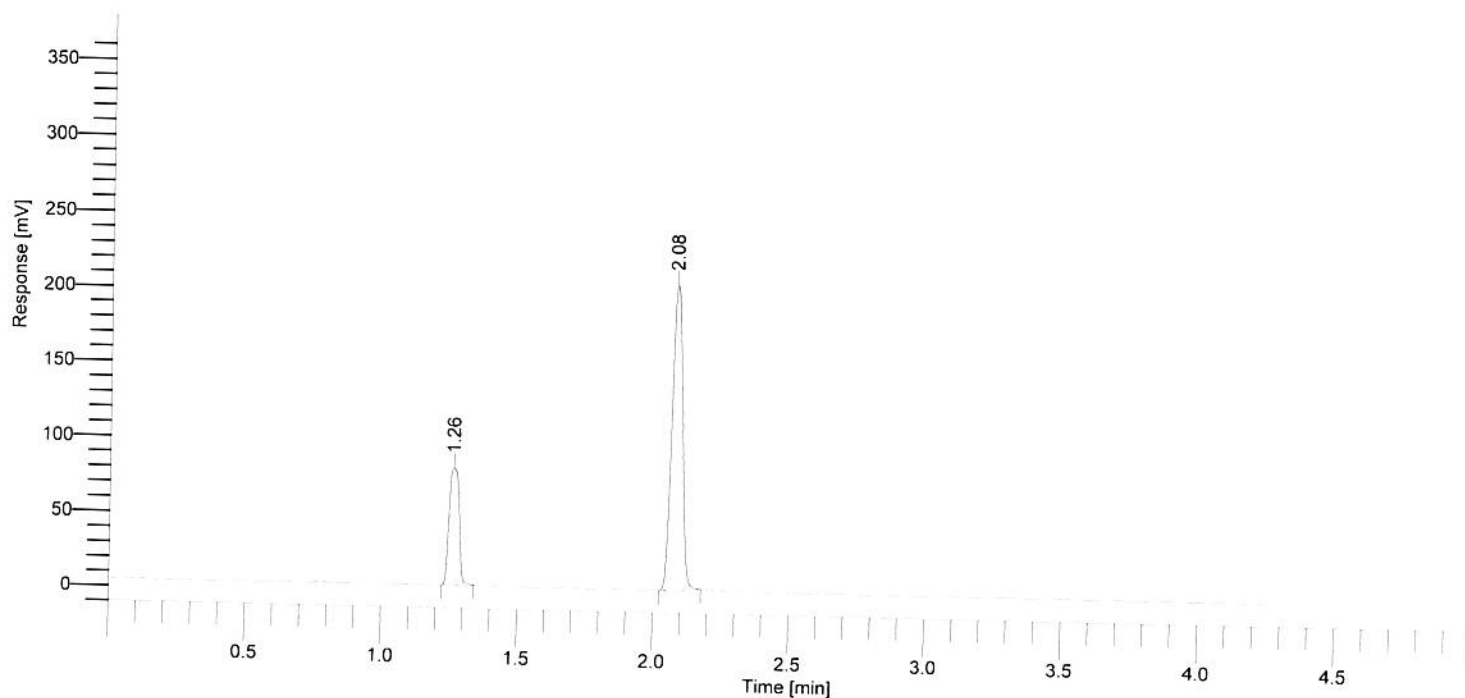

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1194	1.26	195358.5
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	545920.8

46 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33317
Sample Name : GUTH 17310 #1391
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 22

Date : 10/5/2017 4:28:12 PM
Data Acquisition Time : 10/5/2017 4:23:46 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-022.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

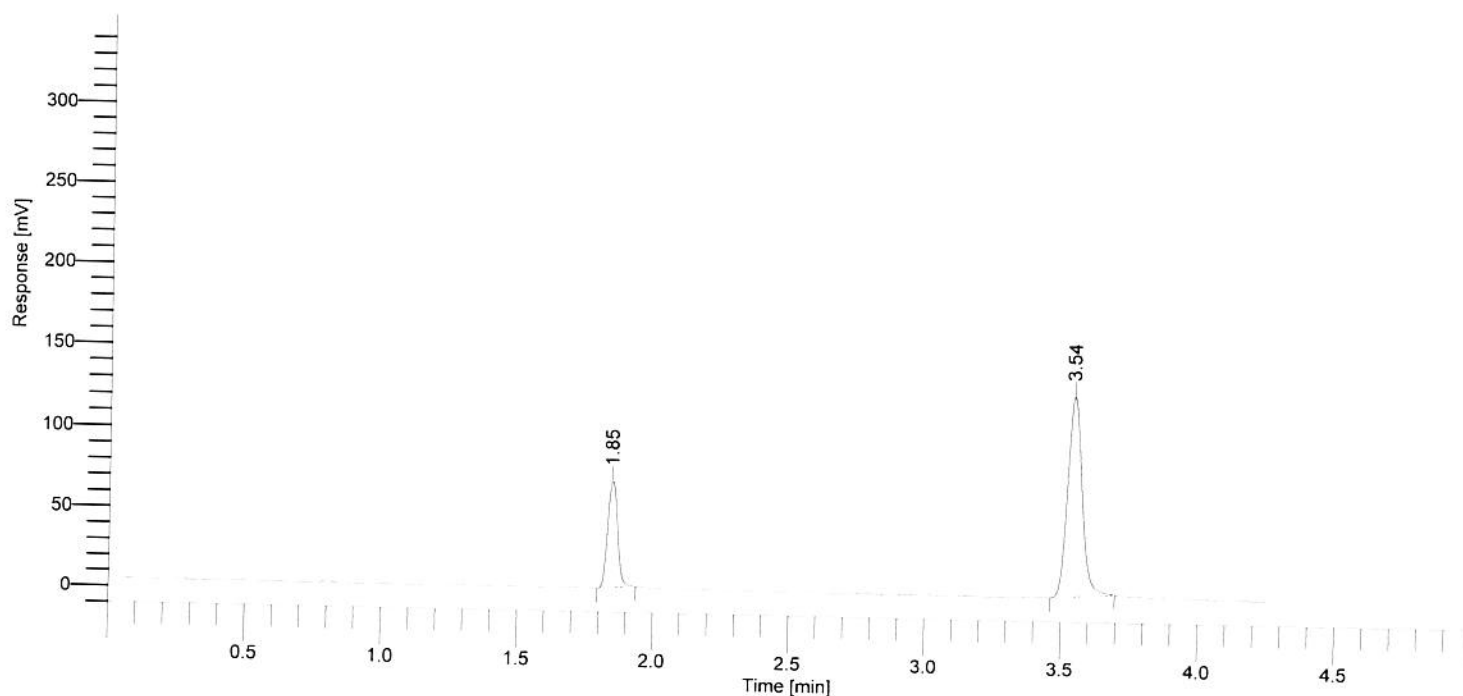

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

47 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33318
Sample Name : GUTH 17310 #1391
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 23

Date : 10/5/2017 4:33:10 PM
Data Acquisition Time : 10/5/2017 4:28:46 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-023.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

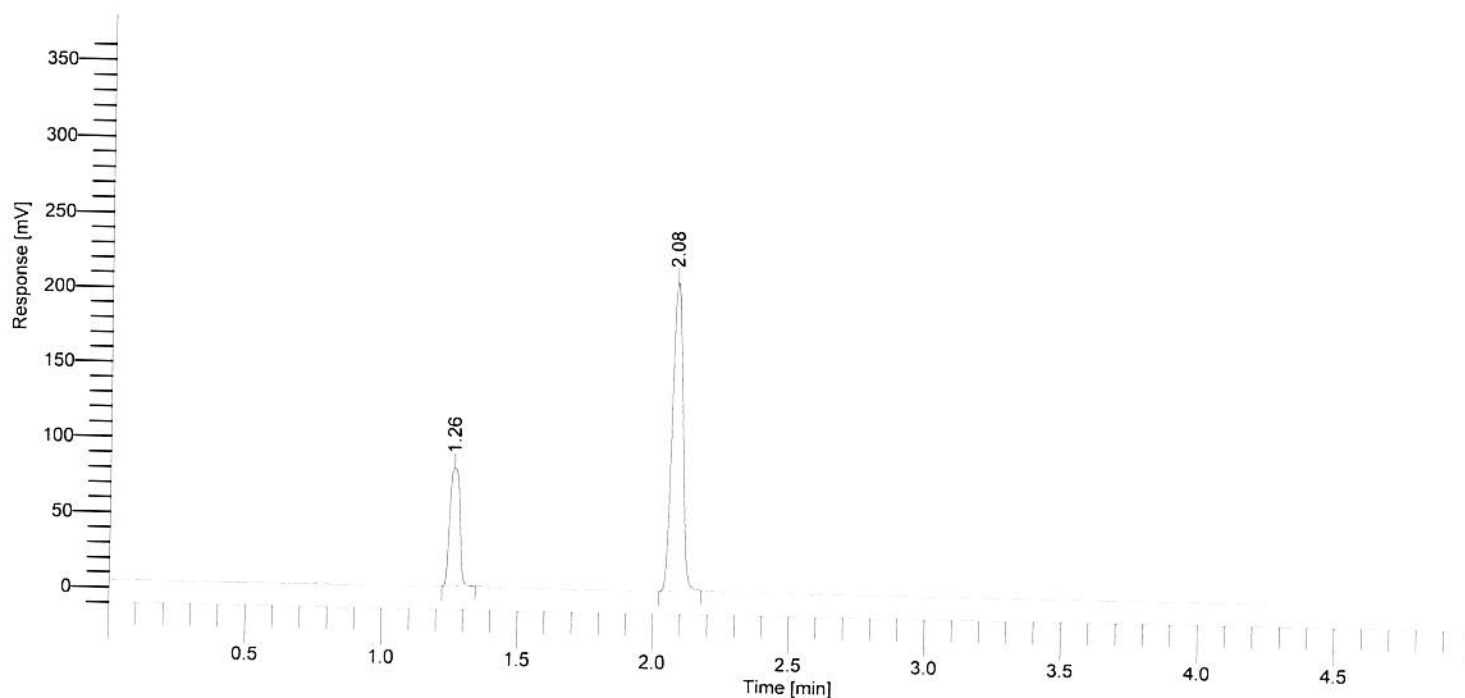

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1195	1.26	197434.8
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.08	551350.4

48

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33319
Sample Name : GUTH 17310 #1391
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 23

Date : 10/5/2017 4:33:12 PM
Data Acquisition Time : 10/5/2017 4:28:46 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-023.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

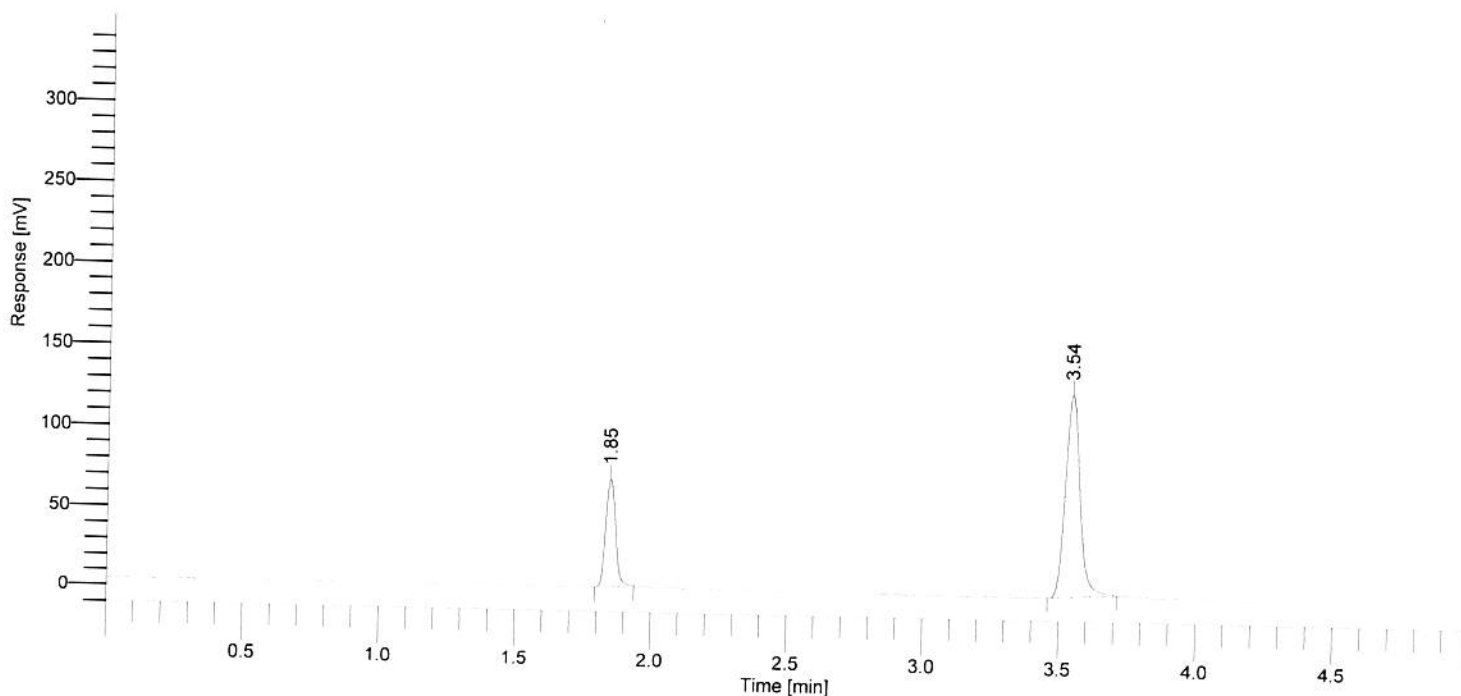

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.85
N-PROPANOL	3.54

49 cm

Software Version : 6.3.2.0646
 Reprocess Number : d-2ua3120y7b: 33320
 Sample Name : GUTH 17310 #1391
 Instrument Name : GCHeadspace2
 Rack/Vial : 0/0
 Sample Amount : 1.000000
 Cycle : 24

Date : 10/5/2017 4:38:11 PM
 Data Acquisition Time : 10/5/2017 4:33:46 PM
 Channel : A
 Operator : manager
 Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-024.rst
 Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

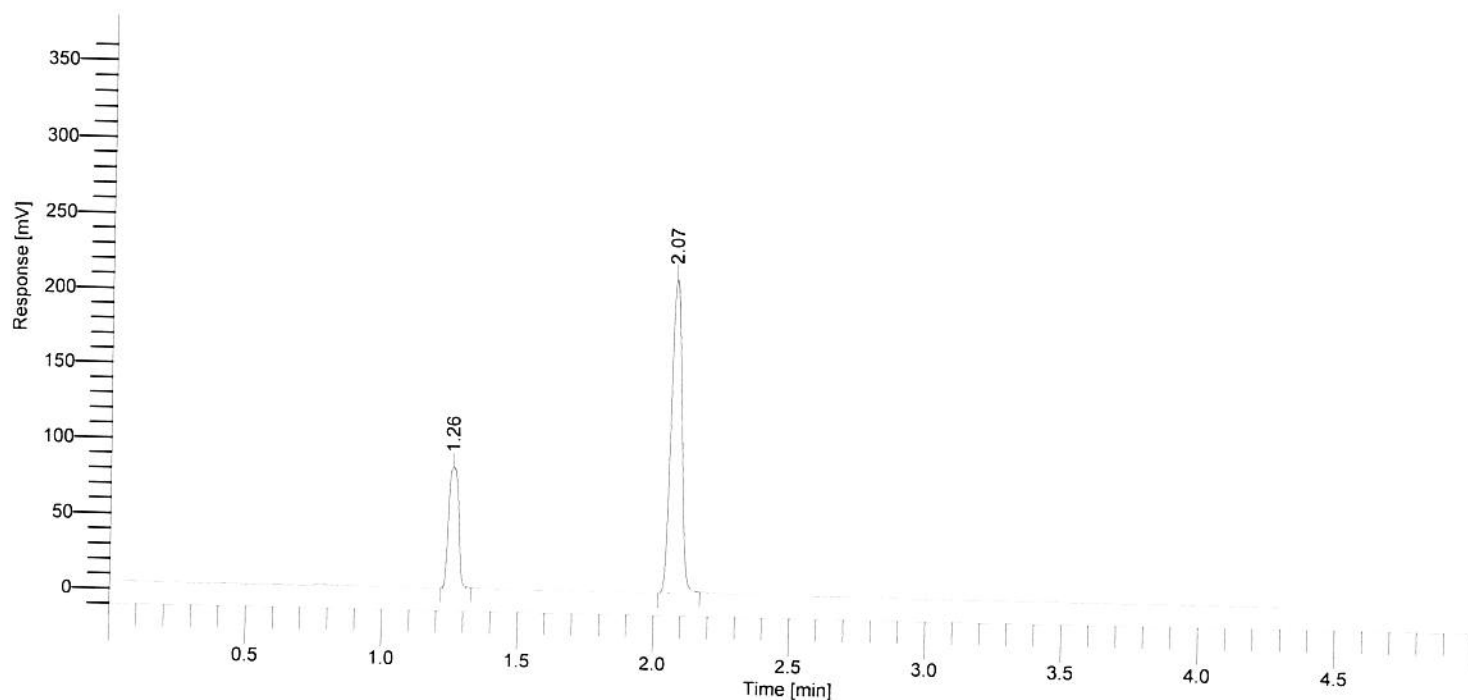

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1188	1.26	199947.7
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.07	561848.2

50 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33321
Sample Name : GUTH 17310 #1391
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 24

Date : 10/5/2017 4:38:13 PM
Data Acquisition Time : 10/5/2017 4:33:46 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-024.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

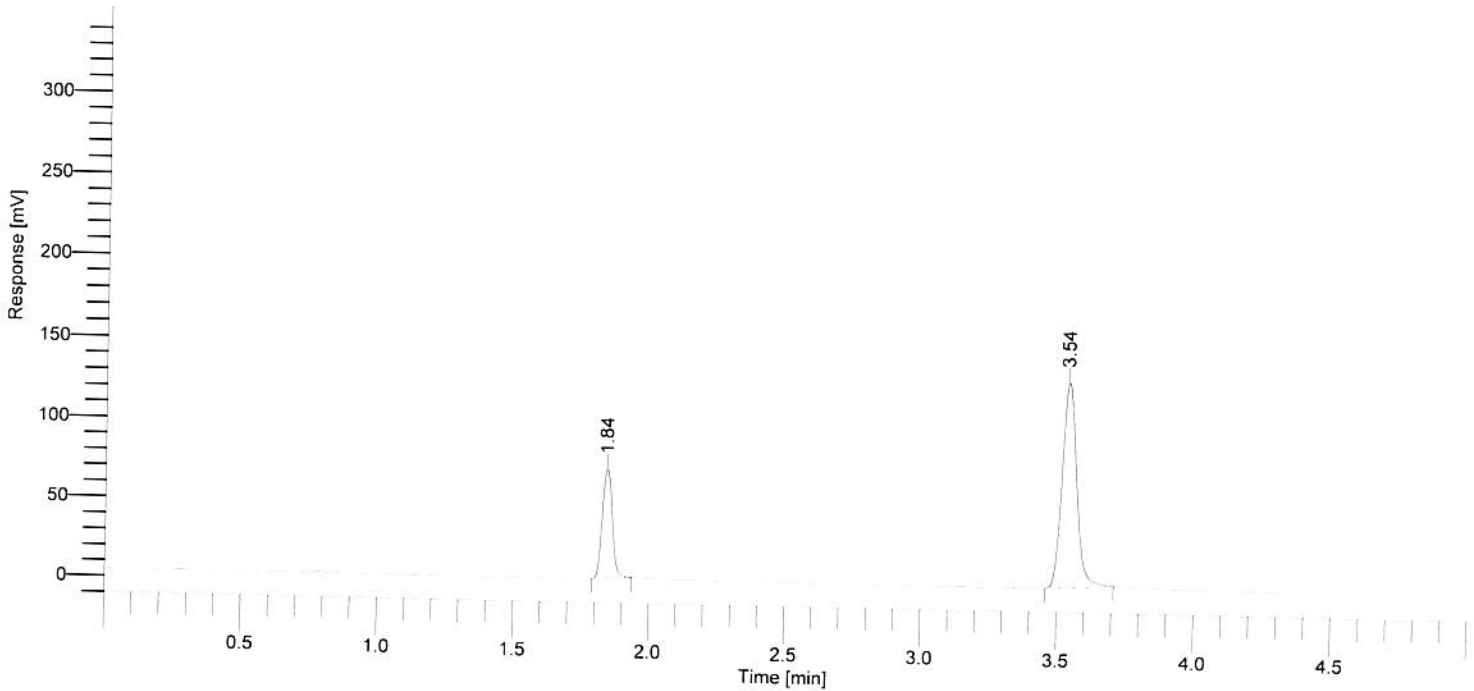

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.84
N-PROPANOL	3.54

SI Cam

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33322
Sample Name : 0.10% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 25

Date : 10/5/2017 4:43:10 PM
Data Acquisition Time : 10/5/2017 4:38:46 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-025.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

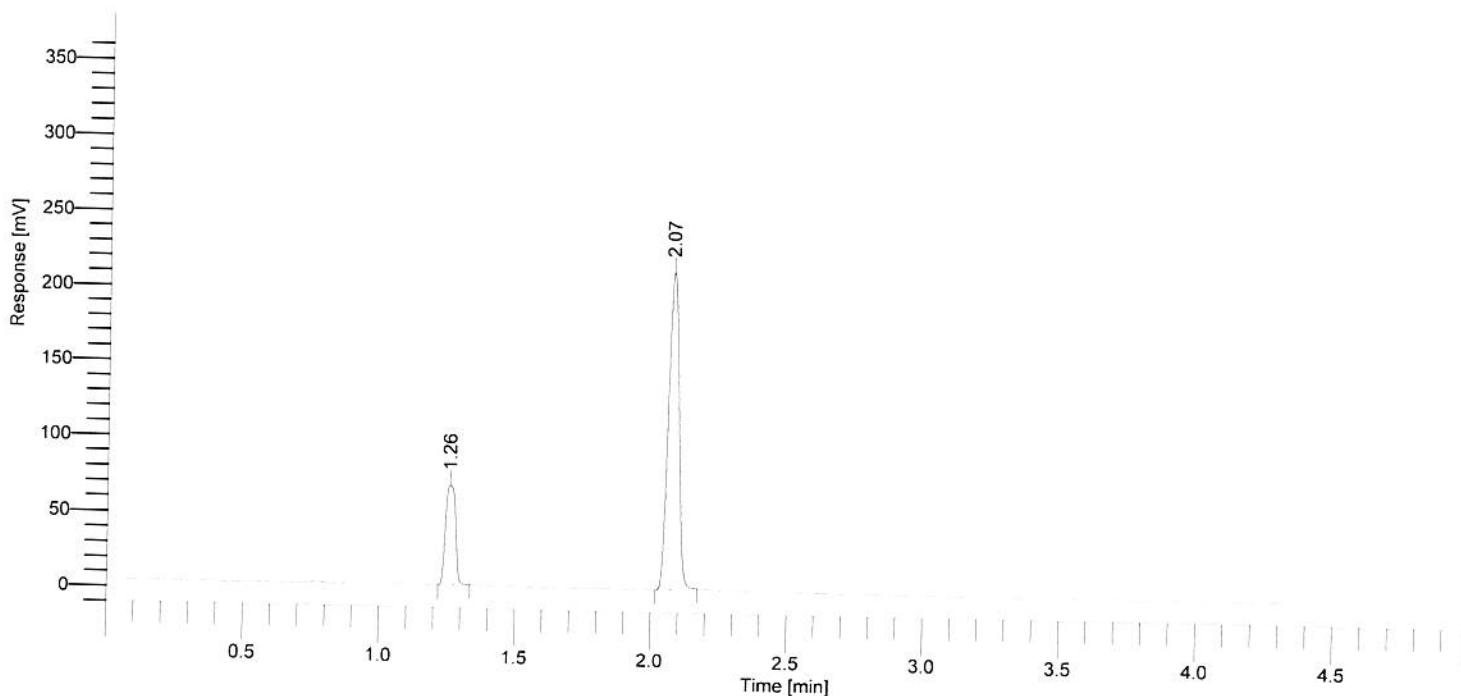

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [$\mu\text{V}\cdot\text{s}$]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.0981	1.26	165818.6
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.07	566306.8

52 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33323
Sample Name : 0.10% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 25

Date : 10/5/2017 4:43:13 PM
Data Acquisition Time : 10/5/2017 4:38:46 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-025.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

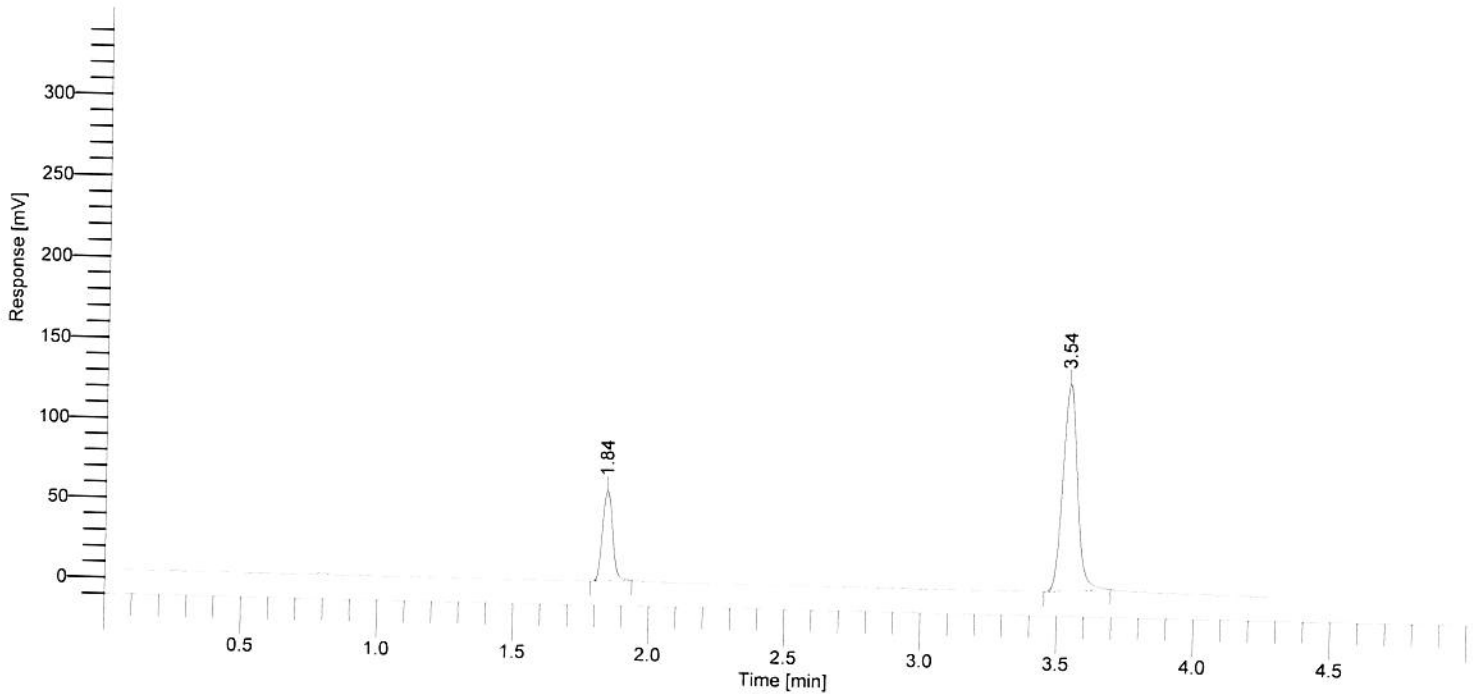

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.84
N-PROPANOL	3.54

53 *Cam*

Software Version : 6.3.2.0646
 Reprocess Number : d-2ua3120y7b: 33324
 Sample Name : 0.20% CONTROL
 Instrument Name : GCHeadspace2
 Rack/Vial : 0/0
 Sample Amount : 1.000000
 Cycle : 26

Date : 10/5/2017 4:48:11 PM
 Data Acquisition Time : 10/5/2017 4:43:46 PM
 Channel : A
 Operator : manager
 Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-026.rst
 Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

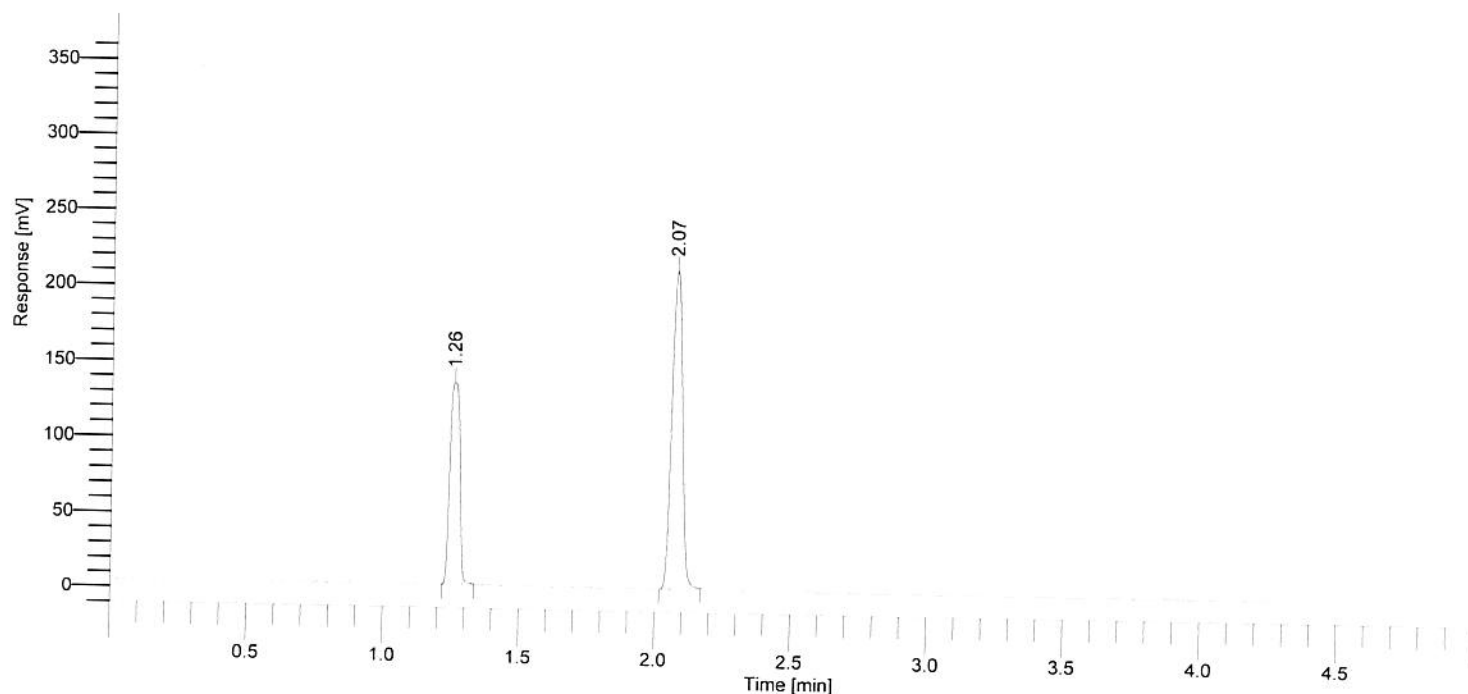

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.1969	1.26	335365.9
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.07	564789.8

54 cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33325
Sample Name : 0.20% CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 26

Date : 10/5/2017 4:48:12 PM
Data Acquisition Time : 10/5/2017 4:43:46 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-026.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

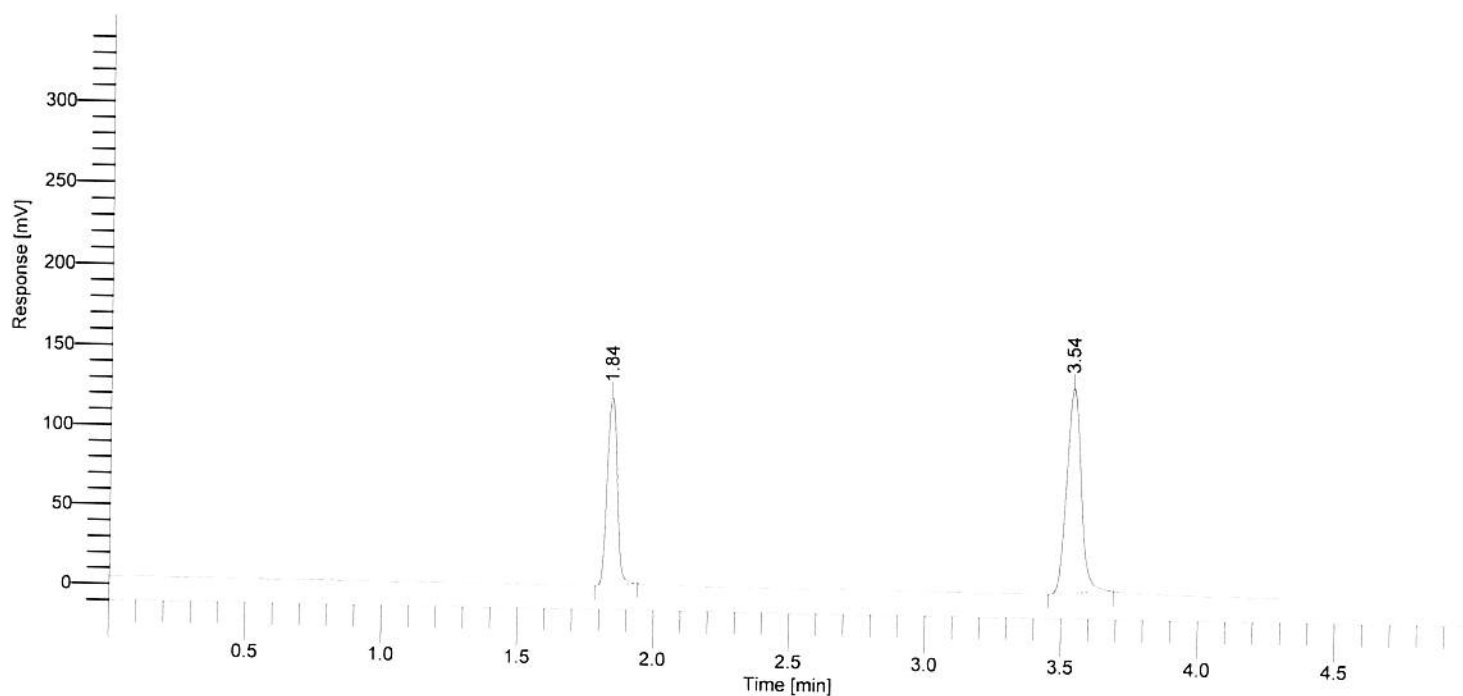

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"


Component Name	Time [min]
ETHANOL	1.84
N-PROPANOL	3.54

SS *[signature]*

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33326
Sample Name : NEG CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 27

Date : 10/5/2017 4:53:11 PM
Data Acquisition Time : 10/5/2017 4:48:46 PM
Channel : A
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\A20171005-1-027.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

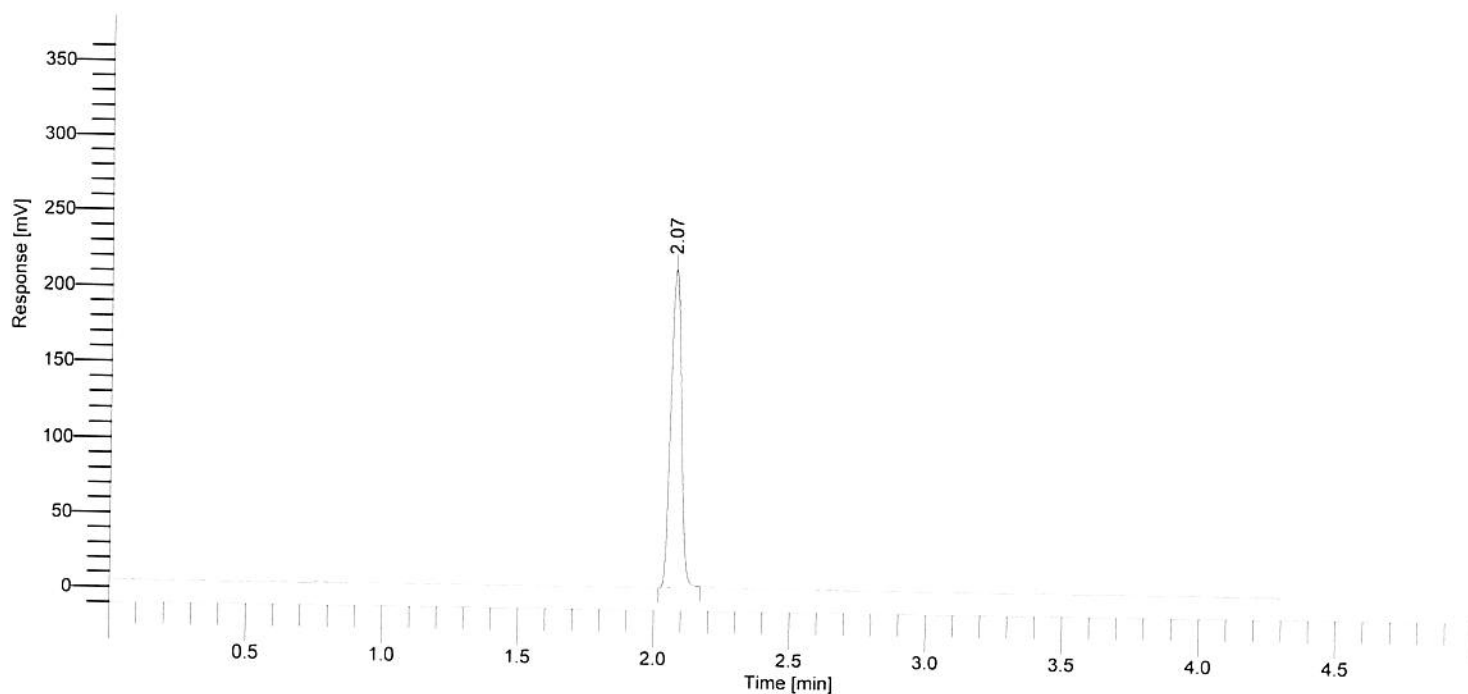

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC1 COLUMN ON FID "A"


Component Name	% by Wt	Time [min]	Area [μ V·s]
METHANOL	0.0000	1.00	0.0
ETHANOL	0.0000	1.26	0.0
ISOPROPANOL	0.0000	1.54	0.0
ACETONE	0.0000	1.90	0.0
N-PROPANOL	-----	2.07	566050.4

SG cm

Software Version : 6.3.2.0646
Reprocess Number : d-2ua3120y7b: 33327
Sample Name : NEG CONTROL
Instrument Name : GCHeadspace2
Rack/Vial : 0/0
Sample Amount : 1.000000
Cycle : 27

Date : 10/5/2017 4:53:12 PM
Data Acquisition Time : 10/5/2017 4:48:46 PM
Channel : B
Operator : manager
Dilution Factor : 1.000000

Result File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Data\B20171005-1-027.rst
Sequence File : C:\PenExe1\TcWS\Ver6.3.2\GCHeadspace2\Sequences\SimSoln.seq

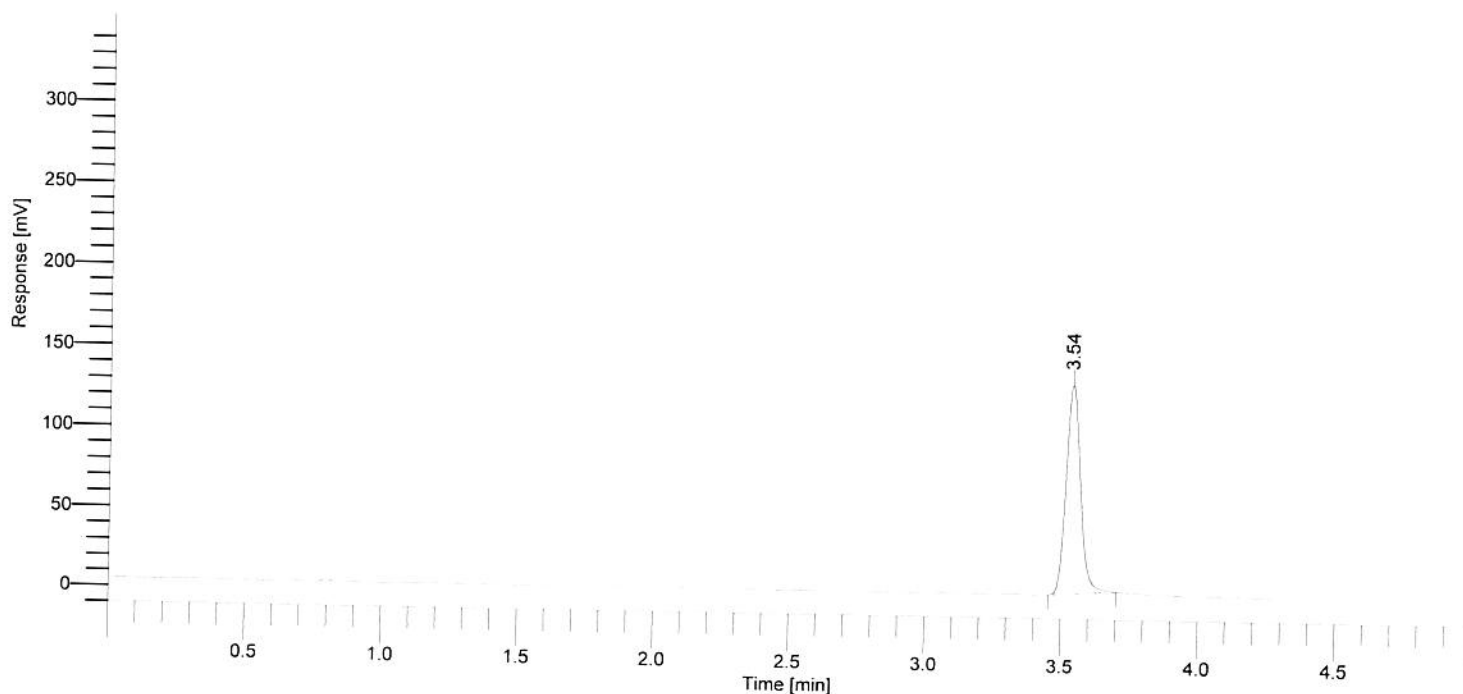

BLOOD ALCOHOL REPORT NYSP

RESULTS FROM BAC2 COLUMN ON FID "B"

Component Name	Time [min]
N-PROPANOL	3.54

S7 Cam